Instant Ear Thermometer

Model 09-340

Please read this instruction manual completely before operating this unit.

INSTRUCTION MANUAL • ENGLISH & ESPAÑOL

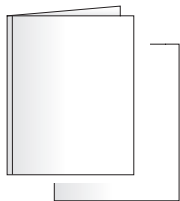
STOP!
PLEASE ENSURE YOU HAVE ALL OF THE
FOLLOWING COMPONENTS BEFORE USING
YOUR INSTANT EAR THERMOMETER.

Thermometer

CR2032 battery

Instruction Manual
Quick Start Guide
FAQ

IF YOU ARE MISSING ANY PARTS, INCLUDING
INSERTS OR INSTRUCTION MANUALS, DO NOT
RETURN TO PLACE OF PURCHASE. CONTACT
CUSTOMER CARE AT 866-326-1313.

TABLE OF CONTENTS

Introduction.....	4
Quick Tips.....	5
Product Features	6
Display Symbols.....	7
Taking a Reading	8
Changing Modes	9
Recalling Readings in Memory.....	9
Replacing the Battery	10
Care and Safety Information.....	11
Cleaning and Storage.....	12
Troubleshooting	13
Device and Label Symbols	14
FCC Statement	15
Product Specifications.....	16
Warranty Information.....	17
Español.....	19-36

Toll-Free Customer Care Help Line: 1-866-326-1313

Monday – Friday 8:30 a.m. – 4:30 p.m. CST

Distributed by:
Veridian Healthcare
1175 Lakeside Drive
Gurnee, IL 60031
www.veridianhealthcare.com

Made in China
#93-1041 05/19
©2019 Veridian Healthcare, LLC

INTRODUCTION

Thank you for purchasing a Veridian Healthcare Instant Ear Thermometer. Please read this instruction manual completely before using your unit, to ensure safe use, accurate results and long life of the thermometer.

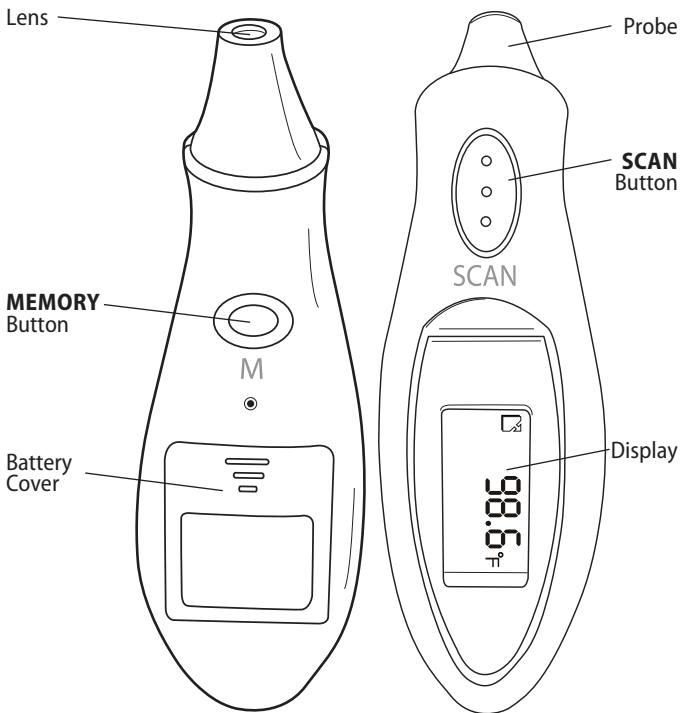
Ear thermometers (also known as tympanic or aural thermometers), are a great alternative to oral, rectal and underarm thermometers. The ear is an excellent site for measuring temperature, as the part of the brain that regulates body temperature shares a blood supply with the eardrum. This thermometer provides results several times faster than the fastest oral thermometer, making it ideal for infants and fidgety children. It can even be used while the patient is sleeping. Unlike readings from oral thermometers, factors such as smoking, drinking, or speaking do not affect the measurements taken by ear thermometers.

The Instant Ear Thermometer takes 512 readings in one second, then converts those data to a single peak temperature, which appears on the display. With proper care, your new thermometer will provide your entire family with many years of safe, accurate and comfortable temperature readings.

QUICK TIPS

- Ear temperature readings are equivalent to oral temperature readings.
- Readings taken while asleep should not be compared directly to readings taken while awake, as body temperature while asleep is typically lower.
- Always take temperature readings from the same ear, as ears may differ in temperature.
- Wait at least 20 minutes after bathing before taking a temperature reading to ensure the most accurate results.
- The best position for infants when taking a reading is lying flat, with the head to the side so the ear being used faces the ceiling.
- For adults and children, a standing position is best. The person taking the reading should stand behind and slightly to the side.
- For certain people, it is best to take three consecutive readings (waiting at least 10 seconds between each), using the highest temperature as the reading. These populations include:
 - Infants less than 100 days old
 - Children under three years with compromised immune systems, or in cases where the monitoring of temperature is important
 - New users of the Instant Ear Thermometer, who are unfamiliar with the technique
- Normal temperature readings can vary from 97.0° F to 99.0° F (36.1°C to 37.2° C). Level of activity, time of day, and other factors may affect body temperature.

PRODUCT FEATURES

DISPLAY SYMBOLS

Measurement in progress

°F / °C

Degrees Fahrenheit
Degrees Celsius

Temperature normal; below
100.4° F (38.0° C); person mode only

Fever alert; temperature
above 100.4° F (38.0° C); person mode only

Low battery

Memory recall

TAKING A READING

1. Clean the lens at the tip of the thermometer, using alcohol and a soft cotton swab. Cleaning the lens before each use ensures the most accurate reading and prevents cross-contamination, see page 12 for complete cleaning and care instructions.

2. Press and release the **SCAN** button. The thermometer will perform a brief self-test, during which all indicators on the display will be shown, Fig. 1, followed by the last reading taken.

Fig. 1

Fig. 2

3. The thermometer will beep, indicating that the unit is ready for measurement. Gently insert the probe into the ear canal until a slight resistance is felt, Fig. 2.

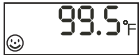
Fig. 3

4. Press and release the **SCAN** button again, while keeping the probe in place. A beep will sound when the measurement is complete.

Fig. 4

5. The temperature reading will be shown on the display.

6. If the temperature measurement is below 100.4° F, the beep indicating the reading is complete will be a single, short beep. A happy face will be displayed next to the reading, Fig. 3.

7. If the reading is 100.4° F or above, the unit will sound a long beep, followed by two short beeps. A sad face will be displayed, Fig. 4.

8. After about one minute after use, the thermometer will automatically beep and shut off.

9. Clean the thermometer lens again after use. Wait at least 10 seconds between consecutive readings, to ensure accuracy.

NOTE: The backlit display will remain lit for approximately 5 seconds each time the **SCAN** button is pressed, after each reading is completed, and each time the **M** (Memory) button is pressed to recall a previous reading.

CHANGING MODES

Your thermometer can display results in either degrees Fahrenheit (°F) or degrees Celsius (°C). To switch between Fahrenheit and Celsius, press and hold the **SCAN** button. While holding the **SCAN** button, press the **M** (Memory) button to toggle between modes. Figs 1, 2.

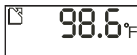
Fig. 1

Fig. 2

RECALLING READINGS IN MEMORY

Your Instant Ear Thermometer automatically stores the last 10 readings taken. To view the readings stored in memory, press the **SCAN** button to turn the thermometer on and wait for the beep. You may then press the **M** (Memory) button to recall the most recent measurement. Press the **M** (Memory) button multiple times to view previous readings (the unit will beep each time you switch readings). The index number to the left of the temperature reading indicates the sequential number of the stored reading you are currently viewing, Fig. 1.

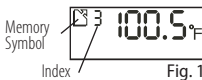
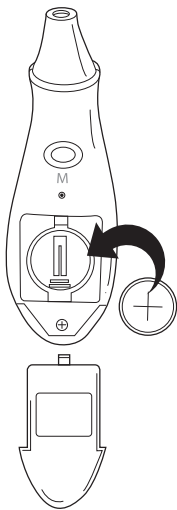
Fig. 1

REPLACING THE BATTERY

Your thermometer came with one lithium battery (type CR2032) installed. When the low battery symbol appears,  the battery needs to be replaced.

Replacing the battery, Fig. 1:

1. Slide the battery cover off
2. Carefully remove the old battery. Do not use a metal-pointed tool.
3. Insert the new battery with the positive (+) side facing up.
4. Replace the battery cover
5. Dispose of the old battery properly at an appropriate facility. Keep the battery away from children, sunlight and heat.

CARE AND SAFETY INFORMATION

Improper care or use of your thermometer may result in injury, damage to the unit or inaccurate readings, which can lead to a misdiagnosed fever. Following these instructions will ensure the thermometer's accuracy and long life.

- Never allow the thermometer to be used without adult supervision.
- Do not use excessive force to insert the probe into the ear; proper use should not cause pain or discomfort.
- The ear should be clean before using this thermometer; earwax in the ear canal may cause an incorrect reading. For proper ear cleaning techniques, consult a healthcare professional.
- Do not use the thermometer if the ear is infected or the ear canal is wet.
- The patient should remain still and avoid talking while the thermometer is in use.
- Clean the lens before and after each use.
- Do not touch the lens with fingers.
- Store the thermometer in the storage case when not in use.
- Do not store the thermometer in direct sunlight, dusty or humid environments, or extreme temperatures. Exact figures for appropriate use and storage can be found in the Product Specifications section of this manual.
- If the patient's temperature is below the ambient air temperature, an inaccurate reading may result.
- Dropping the thermometer or causing other heavy impact may damage it.
- Do not disassemble the thermometer or attempt to repair it (other than replacing the battery).
- If a reading cannot be taken, and the troubleshooting steps do not fix the problem, the thermometer must be replaced.

CLEANING AND STORAGE

The lens should be cleaned before and after each use, to eliminate dirt and fingerprints and ensure the accuracy of your readings. Proper storage will prevent degradation and damage to the unit.

Fig. 1

CLEANING:

- An alcohol swab, or a soft cotton cloth moistened with 70% isopropyl alcohol, should be used to gently cleanse the lens. Do not use water to wash the thermometer.
- The body of the thermometer can be cleaned with a soft, dry cloth.
- Never use a paper towel to clean the probe tip, as it may scratch the lens.
- Do not immerse the thermometer or allow liquid to enter the casing.
- Do not use cleaning agents other than alcohol to clean the thermometer; abrasives, thinners or gasoline will damage the unit.

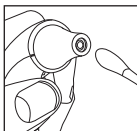
Fig. 2

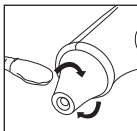
Fig. 3

STORAGE:

- Store the unit in an area away from sunlight, dust and humidity.
- Protect the thermometer from mechanical shock or heavy impact.
- Avoid extended exposure to extreme temperatures. The thermometer should be stored in an environment between -4°F and 122°F (-20°C — 50°C), and with relative humidity of 95% or lower.
- Remove the batteries during extended storage.

TROUBLE SHOOTING

Performance may be degraded should one or more of the following occur:

- Operation outside stated temperature and humidity range
- Storage outside state temperature and humidity range
- Mechanical shock
- Patient temperature is below ambient temperature

The following chart indicates error messages that may appear on your display, why they appear, and how you may solve the problem.

Error Symbol	Problem or Cause	Solution
	Low battery warning.	Replace battery when possible (thermometer will work).
	Low battery.	Replace battery immediately.
 (blank display)	Unit turned off or battery problem.	Press SCAN button to activate unit. Check battery type or polarity (position in the unit). Replace battery.
	Temperature measurement is below typical temperature range.	Gently reposition the probe firmly into the ear canal and try again.
	Temperature measurement is above typical temperature range.	Gently reposition the probe firmly into the ear canal and try again.
	Measurement is below the normal temperature range.	Reposition the probe firmly into the ear canal and try again.

DEVICE & LABEL SYMBOLS

These symbols may appear on your device, instructions or packaging and may vary by make and model.

Symbol	Meaning
	Read This Manual—All included manuals should be read prior to first use.
	Warning—Symbol indicates a warning, prohibition or mandatory action that mitigates a risk that is not necessarily obvious to the device operator.
	Type BF Applied Parts—Indicates that a part of this unit comes in contact with the patient in order to carry out its intended function.
	Environment Protection—dispose of this product properly; consult with your local recycling ordinances for proper recycling and disposal.
	Manufacturer—This symbol indicates the date of manufacture.
	Serial Number

FCC STATEMENT

IMPORTANT INFORMATION REQUIRED BY THE FCC

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference; and (2) this device must accept any interference received, including interference that may cause undesirable operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

Reorient or relocate the receiving antenna.

Increase the separation between the equipment and receiver.

Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.

Consult the dealer or an experienced radio/TV technician for help.

CAUTION: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

SPECIFICATIONS

Name	Instant Ear Thermometer
Model number	09-340
Display system	Liquid Crystal Display with 0.1° F (0.1° C) resolution
Accuracy	+/- 0.4° F, 96.8° F – 102.2° F (+/- 0.2° C, 36.0° C – 39° C) +/- 0.5° F, <96.8° F or >102.2° F (+/- 0.3° C, <36.0° C or >39° C)
Measuring Range	93.2° F – 109.4° F (34.0° C – 43.0° C)
Memory	Automatically stores up to 10 readings.
Tone	One short beep when the unit is turned on Once the measurement is complete: <ul style="list-style-type: none">• One beep sounds for temperatures under 100.4° F (38.0° C)• One long beep followed by two short beeps sound for temperatures above 100.4° F (38.0° C).• One short beep sounds for temperatures <93.2° F (34° C) (Lo) or >109.4° F (43° C) (Hi). Beeps once for each memory recall Beeps once when unit automatically turns off Short rapid continuous beeping indicates battery is low
Operating temperature	60.8° F to 104.0° F (16.0° C to 40.0° C) at relative air humidity (not condensing) or up to 95%.
Storage temperature	-4.0° F to 122.0° F (-20.0° C to 50.0° C)
Automatic power-off	Automatic shut-off after approximately 1 minute of non-use.
Battery	1 each, CR2032 battery
Dimensions	4-1/2" (L) x 1-3/8" (W) x 1" (H) (11.5 cm x 3.5 cm x 2.5 cm)
Weight	With battery, 1.6 ounces (45 grams)
Accessories	Storage case, Instruction Manual, Quick Start Guide, Frequently Asked Questions, one CR 2032 battery
ASTM	E1965–98

Specifications are subject to change without notice.

FIVE-YEAR LIMITED WARRANTY

The warrantor guarantees that its Instant Ear Thermometer will be free from manufacturing defects under normal use for a period of five years for the original purchaser. This warranty covers only normal use, and does not apply to use in any clinical or commercial applications. This warranty does not cover batteries or other power sources that may be provided with or used with the Ear Thermometer. This warranty is voided if the Ear Thermometer product is misused or abused in any manner.

If the Instant Ear Thermometer fails to operate during the time the original purchaser owns it, return it postage prepaid to: **Repair Department, 1175 Lakeside Drive, Gurnee, IL 60031**. The warrantor will repair or replace the defective unit at the option of the warrantor. Repair or replacement of the defective unit, at warrantor's option, is the sole remedy under this warranty.

ANY IMPLIED WARRANTIES WHICH THE PURCHASER MAY HAVE ARE LIMITED IN DURATION TO THE TIME THAT THE ORIGINAL CONSUMER PURCHASER OWNS THE PRODUCT. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

This warranty constitutes the warrantors only responsibility and obligation to repair and/or replace materials or components, or refund the purchase price. The warrantor will not be responsible for any indirect, incidental, special, consequential, or punitive damages or other loss, including, but not limited to, damage to or loss of other property or equipment and personal injuries, whether to purchaser or others. The warrantor shall in no event be liable to the purchaser for any amount in excess of the cost of repair and/or replacement of the unit, or the purchase price of the unit. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

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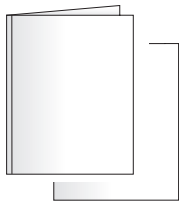
Distributed by:
Veridian Healthcare
1175 Lakeside Drive
Gurnee, IL 60031
www.veridianhealthcare.com
Hecho en China
#93-1041 05/19
©2019 Veridian Healthcare, LLC

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