



SKIN REJUVENATING LASER

Instructions for Use

WELCOME

WELCOME TO LASER SKIN REJUVENATION WITH TRIA BEAUTY. ENJOY PROFESSIONAL ANTI-AGEING RESULTS AT HOME.

With the Tria Skin Rejuvenating Laser you can reverse multiple signs of ageing, restoring your skin's natural radiance. Your new device is a true revolution in personal skincare because it gives you the benefits of professional laser skin rejuvenation through short, simple, daily treatments in the comfort of your own home. The laser light stimulates your body's own natural renewal process, regenerating cells and producing new collagen. The result is more youthful-looking skin with fewer fine lines and wrinkles, a more even colour, and a smoother texture. Proven over a decade of clinical use and favoured by dermatologists, laser skin rejuvenation treats multiple signs of ageing safely and with better results than traditional anti-ageing creams and lotions. You'll have visible results in as few as two weeks.

Before you use your Tria Skin Rejuvenating Laser, please read these *Instructions for Use* and watch the DVD (video also available at www.triabeauty.co.uk/tips) for important safety information and to ensure best results.

TABLE OF CONTENTS


GETTING TO KNOW YOUR TRIA SKIN REJUVENATING LASER	4
BEFORE USING YOUR TRIA SKIN REJUVENATING LASER FOR THE FIRST TIME ...	6
HOW TO USE YOUR TRIA SKIN REJUVENATING LASER	7
HELPFUL TIPS FOR BEST RESULTS	12
HOW OFTEN TO USE YOUR TRIA SKIN REJUVENATING LASER	13
TAKING CARE OF YOUR TRIA SKIN REJUVENATING LASER	14
SPECIFICATIONS	15
TROUBLESHOOTING	16
IMPORTANT SAFETY INFORMATION	19
WARRANTY	22

GETTING TO KNOW YOUR TRIA SKIN REJUVENATING LASER


WHAT'S INCLUDED:

- TRIA SKIN REJUVENATING LASER
- CHARGING CRADLE
- POWER SUPPLY AND CORD
- INSTRUCTIONS FOR USE
- INSTRUCTIONAL DVD

IMPORTANT FEATURES:

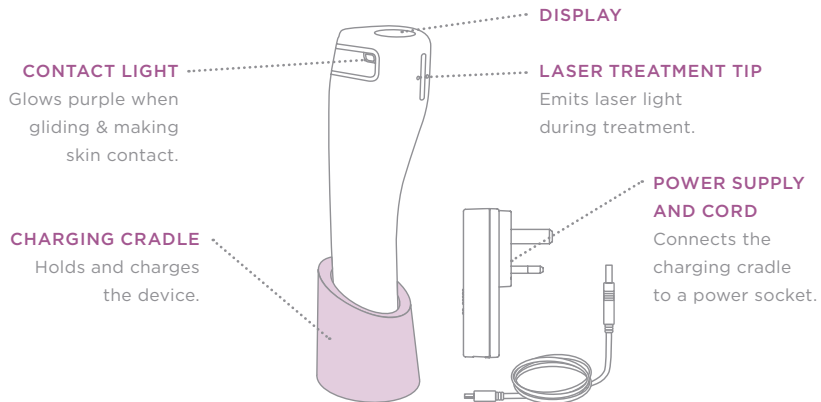

TRIA SKIN REJUVENATING LASER

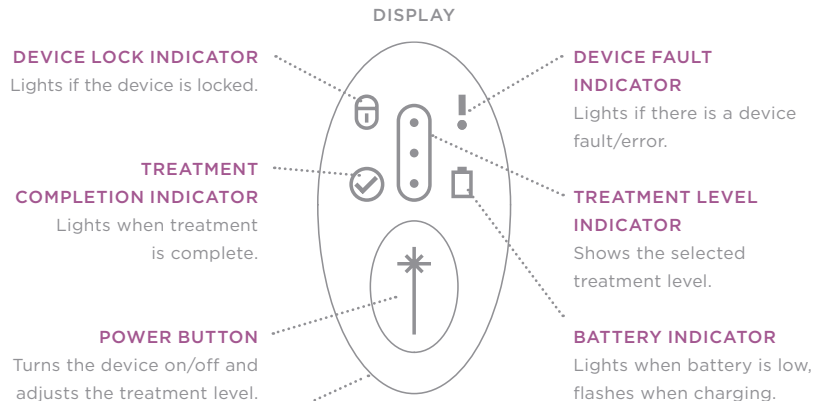

LOCK FEATURE

To lock the device: First, your Tria Skin Rejuvenating Laser must be on. Then press the power button for several seconds until the padlock icon illuminates on the display.


NOTE: When the device is locked it will not power on. The blue padlock icon will light when you press the power button.

To unlock the device: Press and hold the power button until you hear four beeps. Then release your finger from the power button.

BEFORE USING YOUR TRIA SKIN REJUVENATING LASER FOR THE FIRST TIME

CHARGE YOUR TRIA SKIN REJUVENATING LASER

- Insert the cord into the power supply and the charging cradle. Then plug the power supply into a socket.
- Place your Tria Skin Rejuvenating Laser into the charging cradle.
- When your Tria Skin Rejuvenating Laser is charging, the battery indicator will flash. When fully charged (about 2½ hours), the battery indicator will turn off.
- For your convenience, your Tria Skin Rejuvenating Laser may be stored in its charging cradle between treatments. This ensures that it is always fully charged, and ready for treatment.

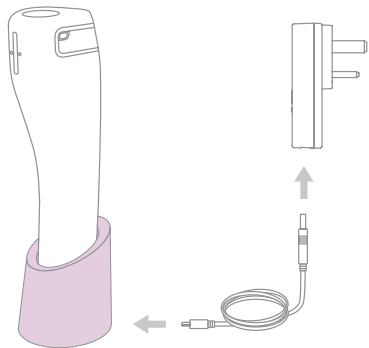

- Please visit www.triabeauty.co.uk/activate or call Tria Beauty Customer Services on 0800 0484 322 to activate.

WATCH THE INSTRUCTIONAL VIDEO

- A one-time activation is required before you use your Tria Skin Rejuvenating Laser for the first time.
- This takes only a few minutes and will also register you for your warranty.

- We've included an Instructional DVD with your Tria Skin Rejuvenating Laser.
- This exclusive DVD gives you helpful information and techniques for using your Tria Skin Rejuvenating Laser.
- Videos are also available at www.triabeauty.co.uk/tips.

ACTIVATE YOUR TRIA SKIN REJUVENATING LASER

HOW TO USE YOUR TRIA SKIN REJUVENATING LASER

THE TRIA SKIN REJUVENATING LASER IS INDICATED FOR THE TREATMENT OF WRINKLES, SKIN DISCOLOURATION AND ROUGH SKIN.


THE INSTRUCTIONS BELOW DESCRIBE HOW TO USE YOUR TRIA SKIN REJUVENATING LASER CORRECTLY.

1. PREPARE YOUR SKIN

- Thoroughly wash and dry your face.
- For best results, use the Tria Rejuvenating Laser Priming Cleanser, which is specially formulated to remove makeup, sunscreen and other products that can block the laser light and make it less effective.


2. SELECT YOUR TREATMENT LEVEL

- Remove your Tria Skin Rejuvenating Laser from its charging cradle and press the power button to turn it on.
- Use the power button to select your desired treatment level (LOW, MEDIUM, or HIGH). For convenience, your Tria Skin Rejuvenating Laser remembers the treatment level you last used.

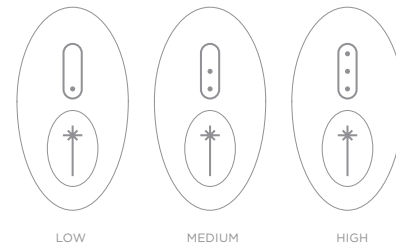

3. TREAT YOUR SKIN

- When treating, divide your face into 4 zones as shown in the picture below to ensure even treatment over the entire face. Treat zone 1 until the device beeps, then move on to zones 2, 3 and 4, in that order.

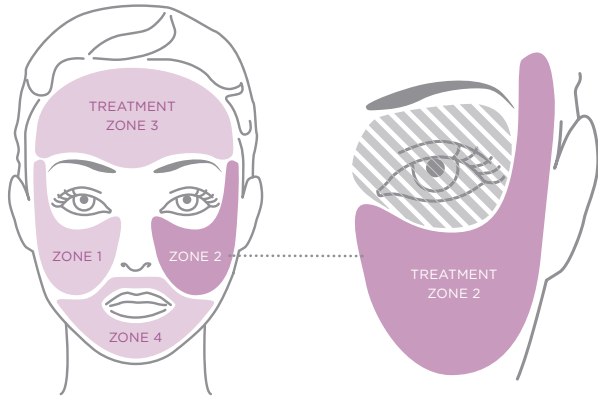

⚠ WARNING: Do not use the Tria Skin Rejuvenating Laser on your eye or eyelid, or within the eye socket. This could cause serious eye injury. Treat around the eyes only where you can feel bone underneath the skin.

3. TREAT YOUR SKIN (CONTINUED)

- To treat each zone, place the entire laser tip flat against your skin and glide it slowly from side-to-side to cover the entire zone as shown in the pictures below.


CORRECT SKIN CONTACT



INCORRECT SKIN CONTACT



CORRECT MOTION

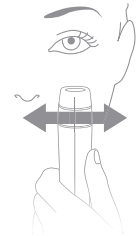

INCORRECT MOTION




Keep the device moving in good skin contact to keep the purple contact light glowing.


3. TREAT YOUR SKIN (CONTINUED)

- When treating your forehead you may find it helpful to hold your device horizontally and move it up and down. Remember to keep the device moving and in good skin contact to keep the purple contact light glowing.

CORRECT
MOTION 1

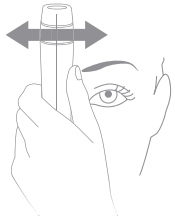

CORRECT
MOTION 2

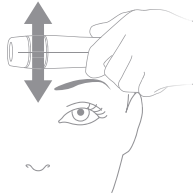

- Treat each zone until your Tria Skin Rejuvenating Laser beeps. This will take about 30 seconds for each zone on LOW, 1 minute for each zone on MEDIUM and 2.5 minutes for each zone on HIGH.
- Your Tria Skin Rejuvenating Laser will beep once after the first zone is treated, twice after the second zone is treated and three times after the third zone is treated.
- After you treat all four zones, the treatment completion indicator (checkmark) will light, a special tone will sound, and the Tria Skin Rejuvenating Laser will automatically turn off.

NOTE: You may feel a slight stinging or warming sensation during treatment. This is normal and generally subsides within 30 minutes. Use of the Tria Rejuvenating Laser Finishing Serum helps reduce these sensations.

NOTE: During treatment or immediately after, you may also develop a temporary redness or slight swelling of your skin. This is normal and most people find that this diminishes over time as their skin adjusts to the treatments. Use of the Tria Rejuvenating Laser Finishing Serum will help soothe your skin.

4. NOURISH YOUR SKIN

- For best results, use the Tria Rejuvenating Laser Finishing Serum, which is specially formulated to improve your laser rejuvenation results by:
 - Soothing and hydrating the skin
 - Exfoliating dead cells
 - Brightening skin tone
 - Reducing pore size
 - Reducing the appearance of fine lines and wrinkles

NOTE: Apply a premium quality sunscreen every morning. With daily sunscreen, you'll protect the benefits from your Tria Skin Rejuvenating Laser treatments, maintain results longer, and help prevent new photodamage (sun damage).

HELPFUL TIPS FOR BEST RESULTS

Watch the instructional DVD

Your DVD contains videos demonstrating how best to hold and move your Tria Skin Rejuvenating Laser across your skin. Videos are also available at www.triabeauty.co.uk/tips.

Practice your technique

Watch yourself in a mirror the first few times you treat. You'll see the glow of the purple contact light when you keep the device flat against your skin and continually moving.

Treat in the evening

Treat in the evening before bedtime to allow any redness or swelling to diminish overnight.

Use the MEDIUM or HIGH treatment level

It may take a week or so for your skin to adjust to the laser treatment. During this period you may want to use the LOW treatment level. Once your skin has adjusted, use a higher treatment level for better results.

HOW OFTEN TO USE YOUR TRIA SKIN REJUVENATING LASER

- For the first 8 weeks, use your Tria Skin Rejuvenating Laser once per day. We recommend making it part of your evening skincare routine. Users typically notice an improvement in as few as 2 weeks with the best results occurring a few weeks after the full 8 week course of treatment.
- Stop using your Tria Skin Rejuvenating Laser for the next 4 weeks to allow maximum cellular and collagen regeneration.
- After this 4 week period, you may begin another 8 week treatment cycle.

TAKING CARE OF YOUR TRIA SKIN REJUVENATING LASER

CLEANING

For best results, keep the laser window and device tip clean. Oils, lotions, makeup, or other debris can block the laser light or keep the device from functioning properly.

If the laser window or device tip does become visibly dirty, clean it with dry cotton wool.

If other parts of the device or charging cradle become dirty, clean as needed with antibacterial wipes.

Make sure the charging cradle is not plugged in and do not get water or other liquid inside any part of the device or the charging cradle.

MAINTENANCE

The device does not need regular or special maintenance. If you think your device or charging cradle is not working properly or is damaged, do not use the device. Call Tria Beauty Customer Services for assistance.

STORAGE

Keep the device and charging cradle away from water, dampness, or debris that may get inside the device or the charging cradle.

DISPOSAL

This device contains lithium batteries. Please follow the regulations where you live for recycling or disposing of electrical devices.

TRAVEL

If you travel with your Tria Skin Rejuvenating Laser, pack the device and charging cradle in its original packaging to prevent damage.

The charging cradle will operate properly at all common household voltages (any supply between 100 and 240 VAC at 50-60 Hz). You may need a plug adapter so that the charging cradle can plug into the wall socket.

LIFE OF THE DEVICE

The device is designed to work properly for years of normal use. There are no parts that need regular replacement. The device is working properly if it functions as described in this booklet.

SPECIFICATIONS

CHARGING CRADLE ELECTRICAL INPUT:

100-240 VAC, 0.2 A max

LASER OUTPUT:

Class 1M, 1440 nm, 12 mJ/pulse max

NOTE: This Class 1M laser device emits invisible laser radiation. Do not treat the eyes, eyelids or point directly at the eye. Do not view directly with optical instruments, such as microscopes and telescopes. This device is safe to eyes and skin if used in accordance with these *Instructions for Use*.

TROUBLESHOOTING

Use this chart to help solve any problems that may arise. If you still have problems, please call Tria Beauty Customer Services on 0800 0484 322.

YOUR QUESTION	IT'S POSSIBLE THAT	HOW TO TROUBLESHOOT
The device fault indicator (!) is lit.	Your Tria Skin Rejuvenating Laser isn't working properly.	Turn your Tria Skin Rejuvenating Laser off, then back on. If the indicator remains lit, contact Customer Services.
My Tria Skin Rejuvenating Laser doesn't charge when I place it in the charging cradle.	Your Tria Skin Rejuvenating Laser isn't connected properly.	Check that the device is sitting correctly in the charging cradle. Then check that the cable is properly connected at both ends and fully plugged into the wall outlet. (see page 6)
	The power socket in your home isn't working.	Try a different power socket.
	Your Tria Skin Rejuvenating Laser isn't sitting in the charging cradle properly.	Re-position your Tria Skin Rejuvenating Laser in the charging cradle, making sure it rests against the bottom of the cradle. Check that the battery indicator is flashing.

YOUR QUESTION	IT'S POSSIBLE THAT	HOW TO TROUBLESHOOT
When I am using my Tria Skin Rejuvenating Laser the contact light doesn't glow purple.	Your Tria Skin Rejuvenating Laser isn't charged.	Charge your Tria Skin Rejuvenating Laser by placing it in the charging cradle.
	The laser treatment tip isn't making good contact with your skin.	Make a slow, steady gliding motion with good skin contact. (see page 9)
	Your Tria Skin Rejuvenating Laser isn't working properly.	Charge your Tria Skin Rejuvenating Laser and restart your treatment session. If the problem persists, contact Customer Services.
When I am using my Tria Skin Rejuvenating Laser the contact light illuminates intermittently.	The laser treatment tip isn't making good contact with your skin.	Make a slow, steady gliding motion with good skin contact. (see page 9)
	Dirt or debris is blocking the laser treatment tip.	Clean the laser treatment tip with dry cotton buds. (see page 14)

TROUBLESHOOTING (CONTINUED)

YOUR QUESTION	IT'S POSSIBLE THAT	HOW TO TROUBLESHOOT
My Tria Skin Rejuvenating Laser won't turn off.	You accidentally released the power button too quickly.	Hold down the power button for 1-2 seconds, then release.
My Tria Skin Rejuvenating Laser won't turn on.	Your Tria Skin Rejuvenating Laser is locked.	Unlock your Tria Skin Rejuvenating Laser. (see page 5)
My Tria Skin Rejuvenating Laser won't lock.	Your Tria Skin Rejuvenating Laser isn't turned on.	First, your Tria Skin Rejuvenating Laser must be on. Then press the power button for at least four seconds until the padlock icon illuminates on the display.
My Tria Skin Rejuvenating Laser is locked/the device lock indicator is lit.	You accidentally locked your Tria Skin Rejuvenating Laser.	Unlock your Tria Skin Rejuvenating Laser. (see page 5)

IMPORTANT SAFETY INFORMATION

Read the following important safety information before using the device. If you have any questions contact Tria Beauty Customer Services (see back page for details) or a health care professional.

WARNINGS

- Do not treat your eye or eyelids. Do not aim the laser directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary vision loss or blindness. Treat around the eyes only where you can feel bone beneath the skin.
- Do not treat skin that is infected, burned or cut. Doing this could cause discomfort or irritation. Wait for these conditions to heal before using the device.
- Do not treat tattooed skin or areas with permanent makeup. Doing this could result in a damaged tattoo or increased risk of skin injury.
- Do not treat the skin more than once per day (once every 24 hours). Doing this could result in increased risk of skin injury.
- Do not treat if your skin is significantly red or swollen from a previous treatment. Wait until the redness and swelling have disappeared before treating again. If you have skin redness or swelling that lasts for more than 3 days, stop using the device and contact a health care professional.
- Do not treat body areas other than the face. Using this device on other areas of the body has not been tested and the risks are unknown.
- Do not allow children to play with or use the device. Doing this could result in significant skin or eye injury.
- Do not use if you are under the age of 18. This device has not been tested on persons under the age of 18 and the risks are unknown.
- Do not use if the device or charging cradle is damaged.

Consult a health care professional before using this device if:

- You have a significant medical skin condition, your skin is unusually sensitive to light, or your skin easily develops a rash or an allergic reaction.
- You are taking a medicine, such as isotretinoin or other oral retinoid (e.g. Accutane), that requires you to avoid sunlight or that contains a steroid, or you have a medical condition that increases sensitivity to light.
- You are undergoing professional cosmetic treatments, such as a skin peel, laser resurfacing or having toxins (e.g., Botox) or fillers injected.
- You are pregnant. This device has not been tested on pregnant women and the risks are unknown.

Failure to follow these warnings may result in potential skin injury.

SIDE EFFECTS

When you are treating, you may feel mild discomfort such as stinging, itching or prickling. There can also be a warming or heat sensation similar to that caused by a sunburn. These sensations are normal and generally subside within minutes or hours and lessen with continued use as the skin becomes accustomed to the treatment. Use of the Tria Rejuvenating Laser Finishing Serum can help minimise these side effects.

Other side effects that may occur in normal use include:

- Skin redness. This typically resolves overnight or within 3 days.
- Skin swelling. This generally diminishes overnight.
- Flaking or dryness.
- Slight temporary darkening of the skin.

If you have more significant or persistent side effects, discontinue treatment and contact Tria Beauty Customer Services or a health care professional.

PRECAUTIONS

To avoid potential device damage or personal injury.

- Do not drop the device on hard surfaces.
- Only use the charging cradle supplied with the device.
- Do not use if the device or charging cradle is wet or could fall into water while it is plugged in.
- Do not attempt to open or repair the device or charging cradle.

Failure to follow these precautions could result in device damage or injury, including vision loss or potentially fatal electrical shock.

WARRANTY

If you have questions about your Tria Skin Rejuvenating Laser or you're not completely satisfied, please refer to the information below and contact Tria Beauty Customer Services.

LIMITED WARRANTY

This document gives you additional legal protection for Tria Beauty products. It is in addition to and not in substitution for your statutory rights as a consumer against Tria Beauty relating to faulty or misdescribed goods or services. Tria Beauty's obligations over and above your normal legal rights to replace faulty product are set out in this document. So please read this document carefully and keep a copy safe. The Tria Skin Rejuvenating Laser comes with a one year limited warranty detailed below.

WHAT IS COVERED BY THE LIMITED WARRANTY?

This one year limited warranty covers defects in materials and workmanship in your Tria Skin Rejuvenating Laser if you purchased it from an authorised partner or directly from the manufacturer (Tria Beauty). Your one year limited warranty is valid for you alone, the original purchaser and cannot be transferred.

WHAT IS NOT COVERED BY THIS LIMITED WARRANTY?

This limited warranty is effective only if the Tria Laser is purchased and used in the United Kingdom. This limited warranty does not cover damage, defects or failure caused by abuse or misuse; improper care; impact or insertion of foreign objects; vandalism; improper storage; improper or incorrectly performed maintenance or repair; application of harmful chemicals; Acts of God, fire, or water; abrasives; negligence; failure to follow manufacturer's instructions; loss of the Tria Skin

Rejuvenating Laser; normal wear and tear; alteration; overuse; professional use or commercial use.

HOW LONG DOES THIS LIMITED WARRANTY LAST?

The limited warranty begins on the date of purchase and expires one year later. The warranty period is not extended if your Tria Skin Rejuvenating Laser is replaced by Tria Beauty during the warranty period.

WHAT WILL TRIA BEAUTY DO IF YOUR TRIA SKIN REJUVENATING LASER IS DEFECTIVE?

If you discover a covered defect in your Tria Skin Rejuvenating Laser before the one year warranty period expires, contact our Customer Services department for assistance. If it is indeed defective, we will replace your Tria Skin Rejuvenating Laser. Customer Services will contact you and keep you updated on the progress of return shipments. You may be required to show proof of purchase (for example a

receipt) prior to the replacement of any Tria Skin Rejuvenating Laser covered by the limited warranty.

WARRANTY LIMITATIONS AND DISCLAIMERS

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND IS IN ADDITION TO YOUR STATUTORY RIGHTS. YOUR REMEDIES FOR THE BREACH OF THIS LIMITED WARRANTY ARE LIMITED SOLELY TO THE REPLACEMENT OF YOUR TRIA LASER AS SET FORTH IN THIS WARRANTY STATEMENT.

IN NO EVENT SHALL TRIA BEAUTY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON BREACH OF WARRANTY, BREACH OF CONTRACT, TORT (INCLUDING NEGLIGENCE), OR OTHERWISE, EVEN IF TRIA BEAUTY IS ADVISED OR AWARE THAT SUCH DAMAGES MIGHT ARISE. SUCH DAMAGES INCLUDE, BUT ARE NOT LIMITED TO, LOSS OF USE OF THE

TRIA LASER, COST OF ANY SUBSTITUTE EQUIPMENT, THE CLAIMS OF THIRD PARTIES, OR DAMAGE TO PROPERTY.

NOTHING IN THIS WARRANTY EXCLUDES OR LIMITS OUR LIABILITY FOR DEATH OR PERSONAL INJURY CAUSED BY A DEFECT IN THE TRIA LASER, OR FOR ANY OTHER LIABILITY THAT CANNOT BE EXCLUDED OR LIMITED BY LAW.

PRODUCT CHANGES

Tria Beauty reserves the right to make changes and improvements to its products without incurring any obligation to incorporate such improvements in products previously sold.

DISPUTES

This warranty and any dispute or claim arising out of or in connection with it or its subject matter or formation (including non-contractual disputes or claims) shall be governed by and construed in accordance with the law of England and Wales. You agree that the courts

of England and Wales shall have non-exclusive jurisdiction to settle any dispute or claim that arises out of or in connection with this warranty or its subject matter or formation (including non-contractual disputes or claims).

TIME FOR FILING CLAIMS

Any claim for breach of this limited warranty must be commenced by within twelve (12) months following the original purchase of the Tria Skin Rejuvenating Laser.

ENTIRE AGREEMENT

Unless modified in writing and signed by both parties, these warranty terms, conditions, and limitations are understood to be the complete and exclusive agreement between the parties, superseding all prior agreements and representations, oral or written, and all other communications between the parties relating to the subject matter of these warranty terms. No employee of Tria Beauty or any other party (including without limitation distributors, dealers, and retailers) is authorised to make any warranty

in addition to those made herein. Nothing in these warranty terms excludes or limits your statutory rights as a customer in respect of faulty or misdescribed goods or services.

MISCELLANEOUS TERMS

As used in this contract, “our”, “we” and “us” refer to Tria Beauty. “You” refers to the original purchaser of the Tria Skin Rejuvenating Laser. A waiver of strict compliance with any term herein shall not be deemed a continuing waiver or waive any other terms or conditions. We shall be excused from performance to the extent of any delay or inability to perform arising due to weather events, disasters, acts of terror, governmental action or regulation, labour or materials shortages or disputes, or other circumstances not under our control. If any term herein is illegal, invalid, or unenforceable, then such term shall be ignored, but the remaining terms and conditions shall remain in effect.

Tria Beauty UK Ltd

Chiltern House
45 Station Road
Henley on Thames
Oxon, RG9 1AT, UK

Customer Services

0800 0484 322
info@triabeauty.co.uk
www.triabeauty.co.uk

