SKINERA™

Beautify your world

IPL Permanent Hair Removal Device 100k


H3230

*with Intense Pulsed Light
Technology*

- (EN) Instructions for use
- (DE) Bedienungsanleitung
- (FR) Instructions d'utilisation
- (NL) Gebruiksaanwijzing
- (ES) Instrucciones de uso
- (IT) Manuale operativo

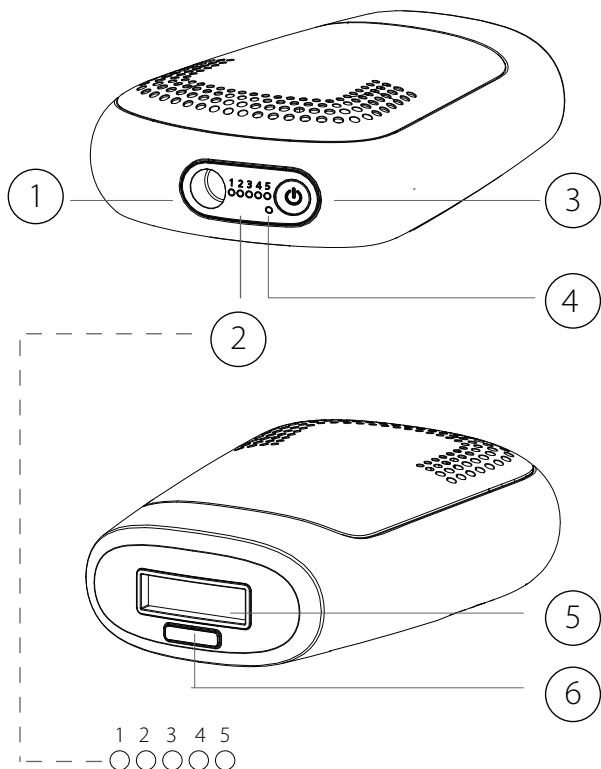IPL Permanent Hair Removal Device 100k

HAIR	SKIN					
		Energy 1-5	Energy 1-5	Energy 1-4	×	×
	Energy 1-5	Energy 1-5	Energy 1-4	Energy 1-3	Energy 1-2	×

Legal Notice

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification.

Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice.

Invention Works B.V., Donk 1B, 2991 LE Barendrecht, The Netherlands,
www.skinera.eu, info@skinera.eu.

Table of contents

EN

1. Get to know your Skinera™ IPL 100k.....	3
2. Safety with Skinera™ IPL device	4
3. Contraindications & warnings	4
4. Do not use your Skinera™ IPL device if:.....	6
5. Consult your physician before use if you have:	6
6. Safety tips.....	7
7. Possible side effects.....	8
8. Long-term hair removal the Skinera™ way	9
9. Get started.....	12
10. After treatment	14
11. Maintenance.....	15
12. Troubleshooting	15
13. Customer service.....	16
14. Warranty	16
15. Disposal	19
16. Specifications	19
17. Signs and symbols	20

1. Get to know your Skinera™ IPL 100k

Skinera™ IPL Permanent Hair Removal 100k is a light-based device for permanent hair reduction designed for home use. It is a powerful electrical device and should be used with special attention to safety. Be sure to read all warnings, contraindications and safety precautions before use, and follow them strictly. We recommend you re-familiarize yourself with the instructions before each treatment.

1.1 Package contents

Your Skinera™ IPL device package includes:

- Hair removal device
- User manual
- Adapter

1.2 Device description

Refer to the diagram at the front of this manual to see the location of each feature of your Skinera™ IPL device.

- | | |
|--|-----------------------|
| 1. Adapter socket | 4. Ready light |
| 2. Energy level indicator lights (1-5) | 5. Treatment surface |
| 3. Control button | 6. Skin colour sensor |

1.3 Intended use

The Skinera™ IPL device is intended to remove unwanted hair and effect long-term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs regrowing after a treatment regime. You may use Skinera™ IPL device to remove unwanted body hair. Ideal body areas include the underarms, bikini line, arms, legs, face, back, shoulders and chest.

Do not use this device for any purpose other than as described in this manual.

2. Safety with Skinera™ IPL device

With Skinera™ IPL safety comes first

Intense Pulsed Light (IPL) technology in the Skinera™ IPL device provides safety with lower energy levels. IPL achieves long-term hair removal results. The low energy reduces the risk of harm or complications, and contributes to your overall safety.

3. Contraindications & warnings

The Skinera™ IPL device is not designed for everyone.

Hair removal by lasers or intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage when used on the face and neck. For more information visit www.skinera.eu.

-  This device must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device. Children should be supervised to ensure that they do not play with the device.
-  Do not use energy levels 2-5 on tanned skin or after recent sun exposure (including artificial tanning machines). Refer to section 6.1.
-  Do not expose treated areas to the sun.
-  Never use around or near the eyes. The device emits infra-red (IR) light and may cause damage.
-  Do not use on areas with injected fills or toxins.

- ❗ Do not treat on tattoos or permanent makeup, dark brown or black spots (such as large freckles, moles, or warts), nipples, genitals, or lips.
- ❗ Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- ❗ Do not treat on an area that has a history of herpes outbreaks.
- ❗ Do not use if you are pregnant or nursing.
- ❗ Do not treat on areas where you may want hair to grow. The results are permanent.
- ❗ Do not treat the same area of skin more than once per hair removal session.
- ❗ Stop use immediately if you see your skin is very red, or if your skin blisters or burns.
- ❗ Use a lower energy level if you experience any discomfort.
- ❗ Remove all jewellery before use.
- ❗ Keep the AC/DC adapter accessible at all times so it can be removed from the socket easily in case of emergency.
- ⊘ Keep your device away from water! Do not use near bathtubs, showers, basins or other vessels containing water. Do not place or store the device where it can fall, be pushed or placed anywhere near water. There is a risk of electrocution. Unplug immediately.

- ! Never leave the device unattended when plugged into an outlet. Keep the AC/DC adapter accessible at all times, so it can be removed from the socket easily in case of emergency.
- ! Do not attempt to open or repair your Skinera™ IPL device. You may be exposed to dangerous electrical components or pulsed light energy. You may cause serious harm to yourself and you will void your warranty.

4. Do not use your Skinera™ IPL device if:

- The plug, attachments or accessories are not recommended by Invention Works B.V.
- The device has a damaged cord or plug. Keep the power cord away from heated surfaces.
- You see or smell smoke when it is in use.
- The device is not working properly, or it appears to be damaged.
- The skin colour sensor is cracked or broken.
- The device casing is cracked or coming apart.

In such cases, contact Skinera™ Customer Service.

5. Consult your physician before use if you have:

- A history of keloidal scar formation.
- Epilepsy.
- A disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- A history of skin cancer or areas of potential skin malignancies.

- Received radiation therapy or chemotherapy treatments within the past 3 months.
- Used physician-prescribed drugs on your skin within the past 6 months.
- Any other condition which, in your physician's opinion, would make it unsafe for you to be treated.

6. Safety tips

This section has information to help you use your device safely.

6.1 Exposure to the sun

- The lowest energy level (1) is designed for safe use even after sun exposure.
- If you use energy levels 2-5, avoid further exposure to the sun for **4 weeks before** and **2 weeks after** your treatment. This is because higher quantities of melanin expose the skin to greater risk when using any method of light-based hair removal. This applies to **all** skin types, even those that don't seem to tan quickly.
- To protect recently treated skin when exposed to sunlight, be sure to thoroughly apply sunscreen SPF 30 or higher, **for 2 weeks after** each hair removal session.

6.2 Exposure of pulse light to the eyes

- Although you can use the device to treat facial hair take extra care to avoid contact with the eyes.
- ! Never use around or near the eyes. The device emits infra-red (IR) light and may cause damage.

6.3 Wounds or infections

- Wait for the affected area to heal before further use.

7. Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important that you understand and accept the risks and complications that may occur with pulsed light hair removal systems designed for home use.

Minor skin discomfort

IPL hair removal is generally well-tolerated. However, some users feel mild discomfort during use. This is best described as a mild stinging sensation on the treated skin areas. The stinging sensation may last during or after treatment.

Anything beyond this minor discomfort is abnormal. This means that either you should stop using the Skinera™ IPL device because you are unable to tolerate the hair removal application, or that the energy level setting is too high.

Skin redness

Your skin may become red immediately after using the device or within 24 hours. Redness generally clears up within 24 hours.

See your doctor if redness does not go away within 2 to 3 days.

Make sure you use the energy level that is comfortable for you and suitable for your skin colour.

Increased sensitivity of the skin

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Pigment changes

The Skinera™ IPL device targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is a risk of temporary hyperpigmentation (increased pigment or brown discolouration) or hypo-pigmentation (whitening) to the surrounding skin. Risk of changes in skin pigmentation is higher for people with darker skin tones or heavily tanned skin. Usually discolouration or changes to skin pigment are temporary, and permanent hyperpigmentation or hypo-pigmentation are rare.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but avoid exposure to the sun.

Blisters or burns

Cool the affected area and apply anti-burn cream. Wait until the area has healed completely before another treatment, and use a lower energy level.

If the skin does not return to normal within 3 days, consult your physician.

8. Long-term hair removal the Skinera™ way

The efficiency of hair removal varies from person to person according to body area, hair colour, and how the Skinera™ IPL device is used. A typical full hair growth cycle may take 18-24 months. During this time, multiple treatment sessions are required to achieve permanent hair removal.

A typical Skinera™ IPL device hair removal plan during a full hair growth cycle includes:

- Treatments 1-4: Plan two weeks apart.
- Treatments 5-7: Plan four weeks apart.
- Treatments 8+: Treat as needed, until desired results are achieved.

8.1 Expected treatment outcomes

Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

Skin Type ¹	Body Hair ²	Light Fluence [J/cm ²]	Part of the Body	Average Treatment Time	Effectiveness of Hair Removal
1-5	Dark blonde to brown	2.5-4	Lower or upper legs	Session 10-12, approx. 4 weeks apart	60%
			Arms		50%
			Bikini line		60%
			Underarms		55%
1-5	Brown to black	2.5-4	Lower or upper legs	Sessions 8-10 approx. 4 weeks apart	70%
			Arms		60%
			Bikini line		70%
			Underarms		65%

¹ Skin types are based on the Fitzpatrick scale of skin colours. Refer to the skin colour chart at the front of this manual.

² Not effective on white body hair.

8.2 Energy levels

The energy level determines the intensity of the light pulse delivered to your skin. As the energy level increases, so do hair removal results, as well as the risk of possible side effects and complications.


The device has 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 *energy level indicator lights* (#2 in the diagram).

When the device is turned on, the energy level is set to level 1 automatically. Only one *energy level indicator light* will be on.

To change the energy level, press the *control button* (#3 in the diagram) repeatedly to increase the energy level from 1 to 5. Once at level 5, press the *control button* again to reduce the energy level from 5 back to 1.

 To determine the right energy level for your treatment area, refer to section "9.1 Before your first treatment".

8.3 Pulsing and gliding

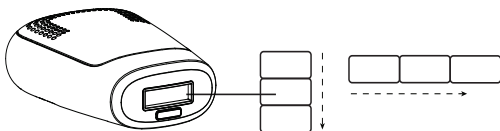Pulsing is the best technique for working at higher energy levels.

Work in rows (either across or down as shown in the image). Start at one end of each row and progress to the other end. Press the *treatment surface* against your skin. The *treatment surface* creates temporary pressure marks on the treated area so you can see areas you have treated. After each pulse, move the treatment surface to the next area. Keep doing this pulsing action for the area you are treating. **Avoid gaps or overlaps between areas.**

 Do not treat the same area of skin more than once per hair removal session.

Gliding is the best technique for working at lower energy levels.

Press the device firmly against your skin and slowly glide over your skin (either across or down as shown in the image). Make sure that the skin is always touching the device. At lower energy levels, the device takes only 1-1.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes.

Glide at a speed so that you avoid gaps or overlaps.

 **Stop use immediately** if your skin blisters or burns.

8.4 The skin colour sensor

Light-based hair removal on darker skin tones may result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypo-pigmentation). Your Skinera™ IPL device has a skin colour sensor that detects the colour of your skin automatically (#6 in the diagram). If your skin tone is too dark for safe use, the device will not operate.

8.5 What to expect when treating with Skinera™

For many people, using the Skinera™ IPL device may be their first experience with a light-based device. During a treatment session it is normal to experience or feel:

- **Noise:** The cooling fan makes noise similar to a hair dryer. This is normal.
 - **A popping sound:** With each flash of IPL light.
 - **A sensation of warmth and tingling:** During each light pulse.
 - **Some redness:** It is not uncommon to see some very mild, pink-like skin colouration during or just after a treatment. This is usually most noticeable around the hairs themselves.
-  **Stop use immediately** if you see your skin is very red.

9. Get started

To use your Skinera™ IPL device, follow these steps:

1. Remove the device and other components from box.
2. Plug the *adapter cord* into the device socket (#1 in the diagram) and plug the *adapter* into an electrical outlet. All energy level indicator lights will turn on, and then go off.

9.1 Before your first treatment

- ! Remove all jewellery before each treatment.

You must decide which energy level is right for your skin and the area you want to treat. The device has 5 levels of energy from 1 (the lowest) to 5 (the highest). To choose the right energy level, test each treatment area separately (refer to the skin colour chart at the front of this manual).

1. Press the *control button* to turn on the device. The level 1 *energy level indicator light* will go on.
2. Place the *treatment surface* (#5 in the diagram) on your skin to deliver one pulse. If you do not experience abnormal discomfort press the control button again. This will increase the energy level. Deliver one pulse on a different spot at energy level 2.
3. Continue testing at increased energy levels. Use the energy levels indicated on the skin colour chart. If you do not experience any adverse effects within one hour, you can start your full treatment at that energy level.

- ! Use a lower energy level if you experience any discomfort.

Repeat this test for each part of the body you wish to treat.

You are now ready to use your Skinera™ IPL device.

1. Before treating an area, the skin should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. **Do not wax, pluck or tweeze the hair out.**
2. Press the *control button*. The device will turn on, the level 1 *energy indicator light* will be on, and a fan sound will start.
3. Approximately 1 second after pressing the *control button*, the *ready light* will turn on (#4 in the diagram). The device is then ready to use at the lowest energy level.
4. **Set the energy level.** Press the *control button* repeatedly to increase the energy level until you have the level you want. The *energy level indicator lights* indicate the energy level setting (refer to section "8.2 Energy levels").

5. Decide if you wish to treat the area by **Gliding** or **Pulsing** (refer to section "8.3 Pulsing and gliding").
6. Apply the *treatment surface* to the skin (#5 in the diagram). Make sure that the skin is spread evenly and smoothly. When the *treatment surface* is in full contact with the skin, the *ready light* will start blinking slowly.
7. The device will first determine the colour of your skin automatically and it will flash a pulse of light onto your skin. You will feel a mild sensation of warmth and tingling, and you may hear a popping sound.

The device will not operate if your skin tone is too dark for safe use.

If you are pulsing, remove the device from your skin. The device will immediately recharge for the next pulse. After a few seconds (depending on the energy level) the ready light will turn on again, and will pulse when pressed against the skin.

8. If you are pulsing, repeat the process. If you are gliding, simply continue.
9. Continue to pulse or glide until you have treated the area fully.

Stop use immediately, if your skin blisters or burns.

10. After treatment

When you have finished your treatment session:

1. Turn the device off by giving the *control button* an extended press. Remember the last energy level setting you used, because it will not be restored when you turn the device on again.
2. Remove the *adapter* from the electrical outlet.

11. Maintenance

It is important to clean your Skinera™ IPL device after each treatment. To clean the device:

-  Unplug the device.
-  To clean the device including the treatment surface, wipe with a clean, dry cloth. Use a cleaning fluid with no more than 70% alcohol, or another non abrasive cleaning fluid. Never use acetone or other abrasive fluids, because they will damage your device.
-  After cleaning, we recommend you store the device in its original box.
-  Never immerse your Skinera™ IPL device or any of its parts in water.

12. Troubleshooting

Do not attempt to open or repair your Skinera™ IPL device. You may expose yourself to dangerous electrical components or pulsed light energy and cause serious harm to yourself. It will also void the warranty.

Problem	Checks
My Skinera™ IPL device does not start.	Make sure the <i>adapter</i> is properly connected to the device and that it is plugged into an electrical outlet on the wall.
There is no light pulse when the device touches my skin.	Check that the device is connected and the ready light is on. Make sure that the treatment surface is evenly and firmly pressed to the skin. It is also possible that your skin colour is too dark.
The skin colour sensor light and the ready light blink.	The skin colour is too dark.
One or more lights blink.	The device has malfunctioned.

For problems or queries, please contact Skinera™ Customer Service.

13. Customer service

See our website for more information: www.skinera.eu. This manual is also available as a PDF download from the website. Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
France	0800-941298	info@skinera.eu
Germany	04473-3270450	info@skinera.eu
Netherlands	0900-0400492	info@skinera.eu
Spain	807001725	info@skinera.eu
UK	0800-3586048	info@skinera.eu
Italy		info@skinera.eu

14. Warranty

What Does This Warranty Cover? Invention Works B.V., its applicable affiliates (hereinafter: Invention Works B.V.) or the authorized dealer from whom you purchased your Skinera™ product (hereinafter: the Device) shall repair or replace the Device, if any material defects caused by Invention Works B.V. are discovered in the Device during the warranty period, whether in the body of the Device or its operation.

Who is entitled to warranty by virtue of this warranty certificate?

Only you, the original buyer of the Device who purchased the Device from Invention Works B.V. or its authorized dealers, are entitled to receive the warranty specified in this warranty certificate. In order to be entitled to the warranty service, you need to provide a proof that you are the original buyer of the Device by presenting an original receipt for purchase of the Device. This warranty is not an international warranty. The warranty will only be honoured in the country in which the authorized dealer or the Invention Works B.V. entity from whom you purchased the Device, are located.

What period does this warranty certificate cover? The warranty covers repair requests received within one year (1) from the Device's original purchase date. Repair requests after one year (1) from the date of purchase are not covered by this certificate of warranty. Repair requests regarding a Device that was repaired or replaced must be submitted no later than: (a) one year (1) from the original purchase date of the Device or (b) ninety (90) days from the date that the repaired or replaced Device was returned to you, the later to occur. If a sticker is attached to your Device stating a warranty period longer than that indicated above, the warranty period mentioned on the sticker will apply, and the Device will be repaired or replaced as part of the warranty until such period expires.

If you purchased the Device in an EU country, you have an additional year (according to the EU-Directive 1999/44) to claim for deficiencies or malfunctioning of your Device, which were present at the time of purchase, but were not apparent then, as determined by Invention Works B.V. or the authorized dealer. In order to make this extended warranty effective you must contact the authorized dealer from whom you purchased the Device or the local service center and notify them about the production fault immediately, and no later than 60 days after detecting the fault.

How to get service for the warranty? For further instructions how to obtain warranty service and for contact details of the service centres close to you, visit www.skinera.eu, or email info@skinera.eu.

When asking for customer services you may be asked to fill in a form for returning the product and you shall receive a reference number. In case of a defect that is not covered by the warranty for whatever reason, Invention Works B.V. can charge you for the shipping and repair costs if you agreed to repair the Device at your expense.

What is not covered by this certificate of warranty? This warranty does not cover any damage resulting from an accident, misuse (including using the Device other than as described in the use instructions accompanying the Device), defective use, lack of reasonable care, irregular and unreasonable use including without limitation commercial or professional use, adding parts not supplied together with the Device, losing parts, using unauthorized spare parts or exposing the Device to any electrical current other than the indicated current.

The warranty will be cancelled if the Device was treated by anyone not authorized by Invention Works B.V. and Invention Works B.V. will not be liable for any damage resulting from repairs performed by anyone not authorized by Invention Works B.V. The warranty does not include replacing used non-durable parts.

WARRANTY DISCLAIMER THE EXPRESS WARRANTY SET FORTH IN THIS DOCUMENT IS IN LIEU OF ALL OTHER WARRANTIES INCLUDING THE IMPLIED WARRANTIES OF MERCHANTABILITY, AND FITNESS FOR A PARTICULAR PURPOSE. THIS WARRANTY IS THE SOLE AND EXCLUSIVE REMEDY FOR ANY DEFECT OF THE DEVICE.

LIMITATIONS IN NO EVENT SHALL INVENTION WORKS B.V., ITS AFFILIATES OR AUTHORIZED DEALERS, BE LIABLE FOR ANY INDIRECT, SPECIAL, PUNITIVE, EXEMPLARY, INCIDENTAL OR CONSEQUENTIAL DAMAGES (INCLUDING, WITHOUT LIMITATION, DAMAGES FOR LOSS OF USE OR CHARGES FOR SUBSTITUTE EQUIPMENT), REGARDLESS OF THE LEGAL THEORY ON WHICH THE CLAIM IS BASED EVEN IF THEY HAVE BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES. IN THE EVENT THAT THE REMEDIES SET FORTH IN THIS WARRANTY ARE HELD TO FAIL OF THEIR ESSENTIAL PURPOSE OR IF INVENTION WORKS B.V., ITS AFFILIATES OR AUTHORIZED DEALERS ARE HELD TO BE LIABLE FOR DAMAGES ON THE BASIS OF ANY LEGAL THEORY, IN NO EVENT SHALL INVENTION WORKS B.V., ITS AFFILIATES OR AUTHORIZED DEALERS, LIABILITY FOR DAMAGES EXCEED THE WHOLESALE PRICE OF A SINGLE UNIT OF THE DEVICE AS OF THE DATE THAT THE CAUSE OF ACTION ACCRUES. YOU ASSUME ALL RISK AND LIABILITY FOR LOSS, DAMAGE OR INJURY TO YOU AND YOUR PROPERTY AND TO OTHERS AND THEIR PROPERTY ARISING OUT OF THE USE OR MISUSE OF, OR INABILITY TO USE, THE DEVICE EXCEPT TO THE EXTENT CAUSED DIRECTLY BY THE NEGLIGENCE OF INVENTION WORKS B.V.

THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS WHICH VARY FROM COUNTRY TO COUNTRY. SOME COUNTRIES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES OR IMPLIED WARRANTIES, SO THE ABOVE DISCLAIMERS AND LIMITATIONS MAY NOT APPLY TO YOU.

15. Disposal

Help protect the environment. Do not dispose with household waste and be sure to follow your country's laws for the disposal of electrical items.

EN

16. Specifications

Model no.	H3230
Spot size	0.9cm x 3cm [2.7cm ²]
Technology	Intense Pulsed Light (IPL)
Maximum energy level	Max 4J/cm ²
Wavelength	475-1200nm
Operation & safety	Skin Colour Sensor ensures use only on appropriate skin types
Package size	180x180x100 mm
System weight	200 gm
Transport & storage between uses and storage condition	Temperature: -40°C to +70°C Relative humidity : 10% to 90% Atmospheric pressure: 500 to 1060 hPa
Operating conditions	Temperature: 10°C to 35°C Relative humidity: 30% to 75% Atmospheric pressure: 700 to 1060 hPa
Time from minimum storage temperature between uses to minimum operating temperature with ambient temperature of 20°C degrees	30 minutes
Time from maximum storage temperature between uses to maximum operating temperature with ambient temperature of 20°C degrees	30 minutes

Mode of operation	Continuous
Adapter models and rating	KSA24A1200150Hy PSM18x-120L6HS For x and y- refer to the country / region
Input Output	AC - 100-240 V, 50-60 Hz, 0.5 A DC - 12 V, 1.5 A
This device was tested and complies with:	Safety: IEC 60335-2-23 EMC: EN 55014 Safety of non laser source: IEC 62471
Expected service life of the equipment	5 years

17. Signs and symbols

	Products marked with this symbol fulfil the requirements of the EU directives.
	Class II. The adapter is protected against electric shock by means of double insulation.
IP30	The device is protected against the ingress of objects larger than 2.5mm diameter, such as fingers, to protect against electric shock.
	Do not use in wet environments.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.

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