

# Metallographic preparation of powder metallurgy parts

Next to casting, mechanical forming and machining, powder metallurgy (P/M) technology is an important method of manufacturing metal parts. Undesirable characteristics of ingot-based metals can be greatly reduced, and desired properties of metals which would normally not alloy easily can be achieved by combining different metal powders or mixtures of metal and non metal powders.

The process of making powders, compacting them into useful shapes and then sintering them is costly, but the finished parts have some specific advantages over wrought or cast parts.

The main advantages are:




- the possibility to make fine grained homogenous structures
- the ability to form complicated shapes with close dimensional tolerances
- and the ability to produce parts with a superior surface finish.

Costly machining processes are thus reduced or eliminated and consequently there is less scrap loss compared to other forming methods. It is therefore most economical to use powder metal-

lurgy for the high volume production of small, intricately shaped, and/or very precise parts such as gears and links.

In addition, the process offers the potential to produce a wide variety of alloys with different material properties such as high temperature toughness and hardness. High speed cutting tool bits from sintered tungsten carbide powder are an example of the variety of different properties which can be achieved with the powder metallurgical process.

As the density of the compacted and sintered part influences its key properties of strength, ductility and hardness, a specific porosity is critical. For process control, metallography is used to check porosity, non-metallic inclusions and cross-contamination. In research and failure analysis, metallography is a major tool used to develop new products and improve manufacturing processes. In addition to chemical analysis, quality control also includes physical methods for checking density, dimensional changes, flow rate etc.

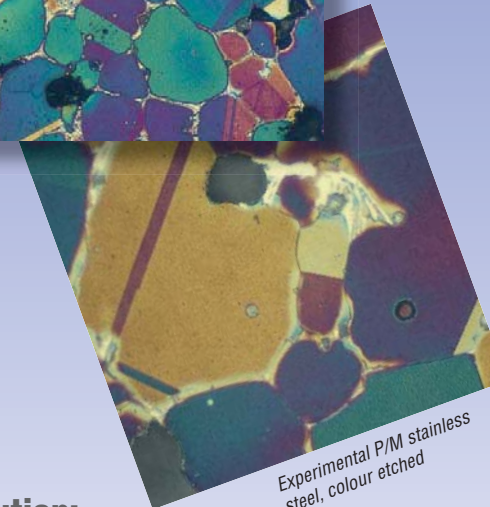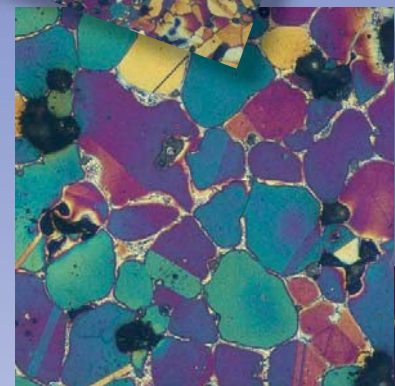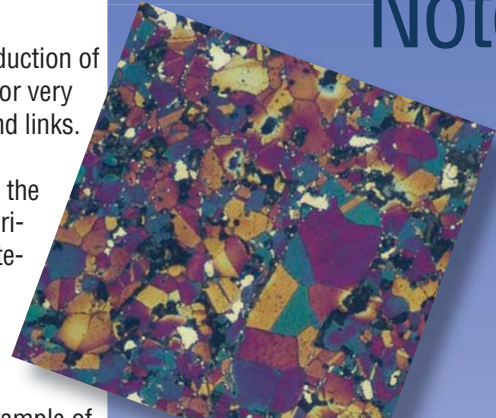

## Difficulties during metallographic preparation

### Grinding and polishing:

To reveal the structure with the correct and representative porosity.



Fig. 1: Porosity of a P/M steel sample after a 4 min. polish with 3 µm diamond



Fig. 2: Same sample as Fig.1 after 8 min. polish with 3 µm diamond

### Solution:

Sufficiently long polishing




Different powder metallurgy steel components

# Production of powder metallurgy parts



Synchronising hub

## Manufacturing process of iron and steel powders

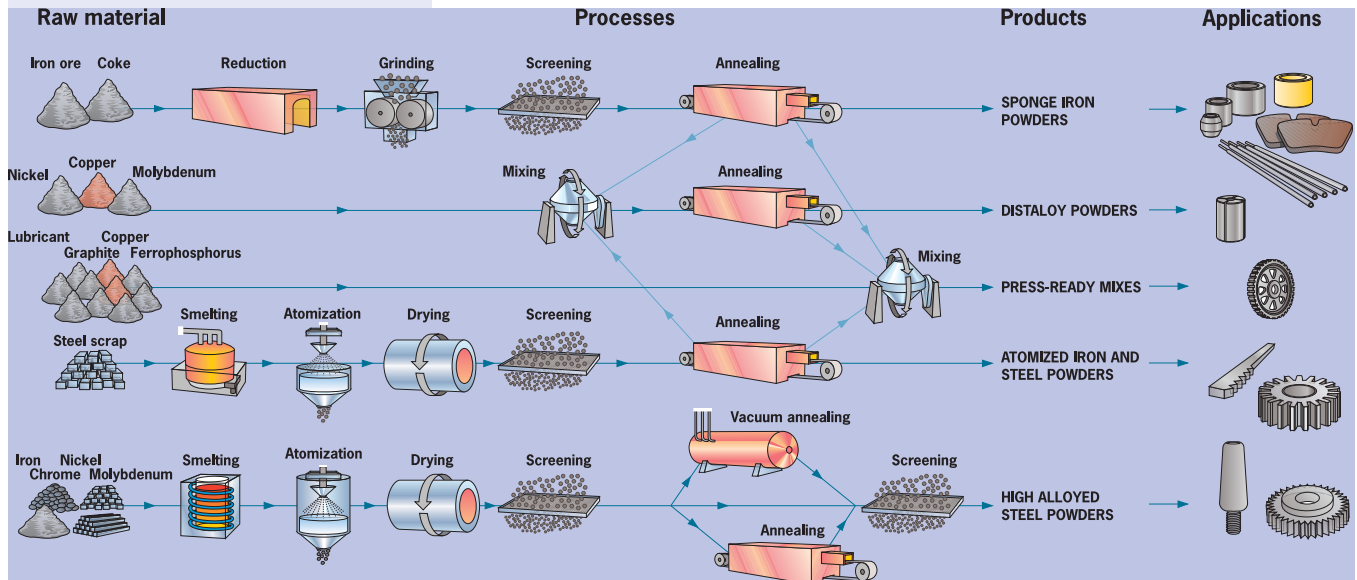

### Production

To achieve the desired structure and near net shape of a powder metal part, stringent process control of the following production steps is required:

- Making the powder
- Mixing the powder with additions such as lubricant, carbon and/or alloying elements
- Compacting the powder in carbide dies
- Sintering at high temperature (1100-1200°C) in a protective atmosphere

Chemical and atomisation methods are the two most common methods for powder production. The chemical method converts metal from ore oxides directly to metal powder at a temperature below the melting point. For example, iron powder is made through direct reduction from iron ore into sponge iron. The sponge iron is then mechanically crushed to powder, which is further refined through annealing under reducing atmosphere to

produce pure iron powder. It is used for alloying and low density applications, such as bearings.

In the atomizing process molten metal of the desired alloy flows through a nozzle and is struck by high pressure water or gas jet. Small droplets are formed which solidify into particles. Atomized powders result in higher densities than mechanically crushed powders, therefore all steel powders are produced by atomizing.

Copper powder is made by atomizing or electrolytic methods. Tungsten carbide powder is produced by adding controlled amounts of carbon to the tungsten powder and carburising it at 1400 - 2650°C


Powder production and mixing is a highly specialized and complex process which produces custom made powder mixes designed to satisfy the needs of a specific application.

A good powder mix not only has the ability to produce the required properties of a specific alloy, but also needs to facilitate handling, compacting and sintering.


For instance, the easy flow of powder and its capability to mix evenly with other powders is important for an even powder distribution before pressing, and ensures uniform properties of the finished part.

For the production of components the mixed powders are first compacted under high pressure in a carbide die. At this stage the part has the geometrical features of the finished component, but not its strength and is called the "green" part. In order to develop the mechanical and physical properties of the material, metallurgical bonding has to take place through sintering at high temperature in a sintering furnace. The bonding occurs through diffusion between adjacent particles. To avoid oxidation, which would impair the inter-particle bonding, the sintering process is conducted in a protective atmosphere.

The bonding increases the density, and pressed and sintered powder metal parts generally contain between 5 and 25% residual porosity. Depending on the application some parts may need an additional hot isostatic pressing for better dimensional accuracy, or surface finish or impregnation with oil. Final treatments, such as surface hardening, plating or coating can be applied.

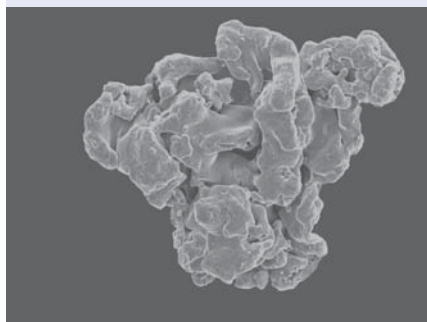

Sponge iron powder, SEM

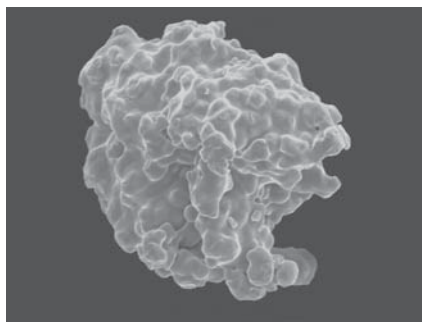

Atomized iron powder, SEM

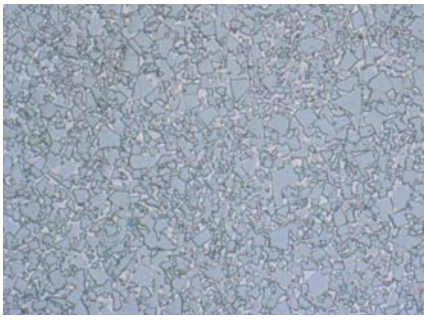


Fig. 3:  
Sintered tungsten carbide  
(WC/Co), etched with  
Murakami's reagent.  
1500 x

The sintering of cemented carbides is carried out in a vacuum sintering process. The carbide powder is mixed with 3-25 w% cobalt and small amounts of titanium and tantalum carbides are added to inhibit grain growth. This mixture is pressed and sintered. At 1280 – 1350°C the liquefaction of cobalt takes place and results in the formation of a eutectic-like phase of WC/Co. Densification begins at lower temperatures and reaches a theoretical 100% shortly after the liquefaction has occurred. During liquefaction the part shrinks up to 40% in volume.

### Applications

Components made by powder metallurgy are mainly used for the following applications:

- Mechanical and structural parts, mainly iron based, but also from copper, brass, bronze and aluminium. The largest user of P/M parts is the automotive industry. Component suppliers make connecting rods, synchronizing hubs, chain sprockets, cams and gears.
- Refractory metals which, due to their high melting points, are difficult to produce by melting and casting.
- Porous material in which controlled porosity serves a specific purpose, for instance self lubricating bearings.
- Composite materials that do not form alloys, for instance copper/tungsten for electrical contacts, cemented carbide cutting tools (Fig. 3), materials for brake linings and clutch facings, diamond cutting tools, or metal matrix composites.
- Special high-duty alloys, such as nickel and cobalt based super alloys for jet engine parts, and high speed tool steels, which have an even distribution of carbides and have isotropic qualities (Figs. 4 and 5).

In addition, different powders and powder mixtures for thermal spray coatings are produced and are also subject to metallographic quality control.

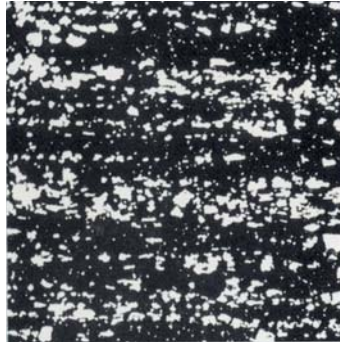


Fig. 4: Carbide distribution in conventionally produced steel



Fig. 5: Carbide distribution in powder metallurgically produced steel



## Difficulties in the metallographic preparation of powder metallurgy parts

The main challenge during the preparation is to show the true porosity after grinding and polishing. Depending on the hardness of the material this can be more or less successful. During the grinding of soft metals, abraded metal is pushed into the pores, which then has to be removed by polishing. Samples from parts in which hard and soft materials are mixed are prone to show pronounced relief. The preparation of green parts needs particular care and patience as they are very fragile.

### Recommendations for the preparation of powder metallurgy parts

#### Cutting

For sectioning a powder metallurgy component of a specific metal or alloy one can select an appropriate cut-off wheel using the recommended charts and guidelines. For mixed materials it is recommended to use a cut-off wheel suitable for cutting the material which constitutes the major part. For sintered carbides a resin bonded diamond cut-off wheel is recommended (e.g. B0D31). Green parts need to be mounted in cold mounting resin (see section "Mounting") before cutting, so that they are not crushed by clamping.

#### Mounting

In order to assure good adhesion of the mounting resin to the sample material it is essential to degrease the sample thoroughly with acetone, toluene or Isopar C\* before mounting (Use proper safety precautions when handling solvents!).

Sintered parts can be hot compression mounted with a resin suitable to the hardness of the sample material, either phenolic resin (MultiFast) or reinforced resins (DuroFast, IsoFast).

Green parts need to be re-impregnated after sectioning under vacuum with a cold mounting epoxy resin (EpoFix, Specifix-40).

Powders can be mounted by mixing a small amount (about 1/2 teaspoon) of powder with a slow curing epoxy resin and pouring it into a mounting cup. During the 8 hours curing process the particles will form a layer by settling at the bottom of the cup.

Hard metal powders can be hot compression mounted by mixing with one measuring spoon of the fine grained mounting resin

\*Available from Exxon

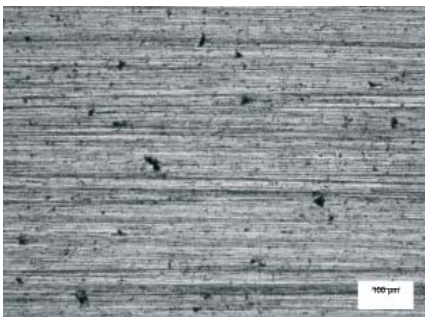


Fig. 6: Surface of P/M steel after fine grinding on MD-Allegro

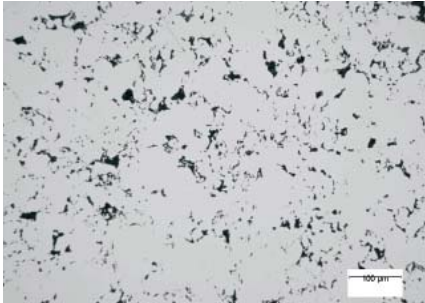


Fig. 7: Same sample as in Fig. 6 showing insufficient polish



Fig. 8: Same sample as in Fig. 7 after longer polish showing correct porosity

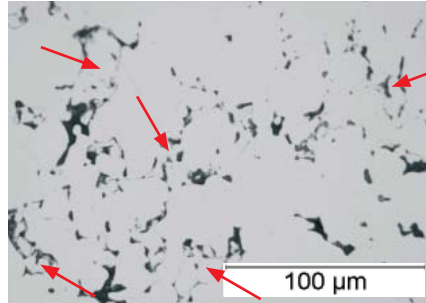


Fig. 9: Larger magnification of surface from Fig. 7, showing metal "lids" covering pores

reveal the correct porosity. For soft metals the final polishing time should not be prolonged unnecessarily, as this leads to rounding of the pore edges. Starting with 500# or 800#, green parts are ground manually on silicon carbide paper to 4000#. The surface is re-impregnated as needed. Polishing can be done on a semi-automatic polishing machine for individual samples.

In order to establish the exact polishing time for specific alloys and parts, it is recommended to check the structure during polishing every two minutes with the microscope, and only proceed to the next polishing step when all residual metal has been removed from the pores. In general, polycrystalline diamond suspension is recommended for polishing

IsoFast. The mixture is then poured into the cylinder of the mounting press and is topped with phenolic resin.

### Grinding and polishing

The grinding and polishing procedure for powder metals follows the same rules which are applied to prepare ingot based samples of the same material.

Plane grinding high volumes of samples of materials >150 HV can be carried out on an aluminium oxide grinding stone or a diamond grinding disc (MD-Piano). Materials <150 HV can be plane ground on silicon carbide paper.

For fine grinding with diamond the fine grinding disc MD-Allegro is suitable for materials >150 HV and MD-Largo for materials <150 HV. This is followed by a thorough diamond polish with 3 µm and a brief final polish with 1 µm or with oxide polishing suspensions. As one of the major goals of preparing a powder metallurgical sample is to show the true porosity, it is important that the diamond polishing step is carried out long enough to achieve this goal (see Figs. 6-9). For large or soft samples it can take up to 10-15 minutes of diamond polishing to remove the metal pushed into the pores during grinding and

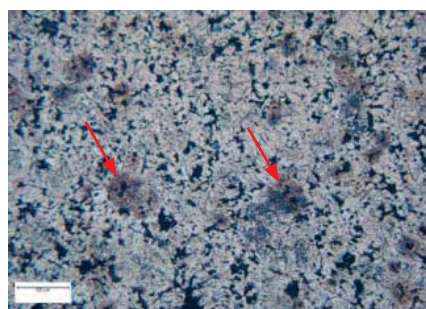


Fig.10: Water stains from cleaning can lead to misinterpretation of structure

powder metals. If the polishing takes excessively long, DiaPro diamond suspension can be used instead.

### Cleaning and drying

After polishing it is essential to clean the sample with a water/detergent mixture to remove remnants of the polishing suspension and lubricant from the pores. The sample is then rinsed with water, followed by a thorough rinse with isopropanol and dried with a

Grinding			
Step	PG	FG	
Surface	SiC paper 320#	MD-Largo	
Suspension		9 µm	
Lubricant	Water	Green	
rpm	300	150	
Force [N]	210	150	
Time	As needed	4 min.	
Polishing			
Step	DP 1	OP	
Surface	MD-Mol	MD-Chem	
Suspension	3 µm	OP-S	
Lubricant	Red		
rpm	150	150	
Force [N]	150	150	
Time	4 min.	1 min.	

Table 1: Preparation method for 6 samples of P/M bronze, mounted, 30 mm dia., clamped into a sample holder, using the semi-automatic TegraSystem.

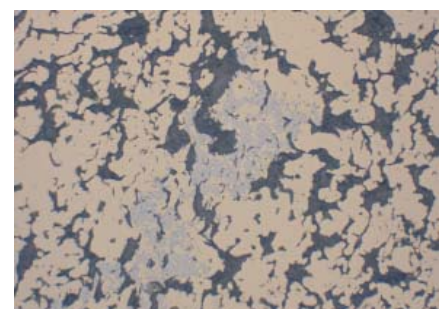


Fig.11: P/M bronze, unetched, containing graphite (grey), and α-δ eutectoid (blue) 500x

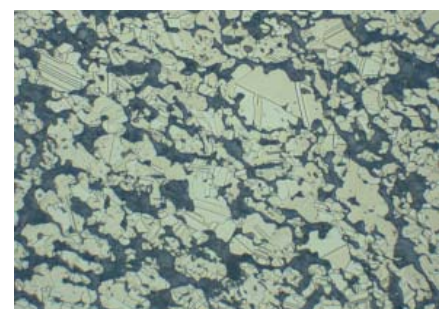


Fig.12: Same as Fig.11, etched with iron-III-chloride, showing grain structure of bronze 500x

Grinding			
Step	PG	FG	
Surface	SiC paper 220#	MD-Allegro	
Suspension		9 µm	
Lubricant	Water	Blue	
rpm	300	150	
Force [N]	210	180	
Time	As needed	5 min.	

Polishing			
Step	DP 1	DP 2	
Surface	MD-Dac*	MD-Nap	
Suspension	3 µm	1 µm	
Lubricant	Blue	Blue	
rpm	150	150	
Force [N]	180	120	
Time	6-8 min.	0.5-1 min.	

Table 2: Preparation method for 6 P/M steel samples, mounted, 30 mm dia., clamped into a sample holder, using the semi-automatic TegraSystem.

\*Alternatively MD-Mol cloth can also be used.

warm stream of air, holding the sample at an angle. Do not blow the stream of air directly from the top onto the sample surface as it forces the liquids out of the pores, which will leave stains on the surface (Fig.10). It is important to use a high quality alcohol for rinsing to keep staining to a minimum.

### Etching

It is recommended to examine the unetched sample first to check the density, shape and size of the pores, oxidation and inclusions, sintering necks and free graphite (see Figs.11 and 12). It is important to know the theoretical density in order to compare it to the porosity.

After this initial inspection it is recommended to etch the sample immediately to avoid drying stains, which form when cleaning and drying liquids are gradually released by the pores. For very dense sintered carbides this

Preparation method for P/M steel

Preparation method for sintered carbides

Grinding			
Step	PG	FG	
Surface	Piano 220#	MD-Allegro	
Suspension		DiaPro Allegro/Largo	
Lubricant	Water		
rpm	300	150	
Force [N]	180	180	
Time	As needed	4 min.	

Polishing			
Step	DP	OP	
Surface	MD-Dac	MD-Chem*	
Suspension	DiaPro Dac	OP-U	
rpm	150	150	
Force [N]	180	120	
Time	4 min.	1 min.	

Table 3: Preparation method for 6 samples of sintered carbides, mounted, 30 mm dia., clamped into a sample holder, using the semi-automatic TegraSystem.

\*Optional step

is not as important as for any powder metal with a certain amount of porosity. To reveal the structure, the common chemical etching solutions for metals and their alloys, that are recommended in the literature, can be used. The following procedure is recommended for etching: wet the surface with isopropanol, immerse sample face up into etchant and slightly agitate. When the appropriate etching time has elapsed, take the sample out of the etchant and rinse either with isopropanol or water depending on the etchant (see below) and dry with a stream of warm air. Interpretation and photographic documentation should take place immediately after drying.

Etching time depends on the alloy and it needs some experience for the correct timing. Etching too short will not give sufficient contrast of the different phases. If the sample is over-etched, it is difficult to distinguish the various phases (see Figs.15-17). When working with unfamiliar material it is recommended to etch shorter rather than longer and to check the result with the microscope first. More etching can be carried out if necessary, but if the sample is overetched, it needs to be repolished. Following are some of the common etchants. Follow standard safety precautions when mixing and working with chemical reagents!

### Copper and copper alloys:

- 100 ml water  
20 ml hydrochloric acid  
5 g Iron-III-chloride  
10-20 sec  
Rinse with water followed by isopropanol
- 100 ml water  
10 g ammonium persulfate. Use fresh only!  
Rinse with water followed by isopropanol

### Steel:

- 1-3% Nital for iron-carbon alloys, iron-carbon-copper alloys and pre-alloyed iron-molybdenum:  
100 ml ethanol  
1-3 ml nitric acid  
10-60 sec depending on carbon content.  
Rinse with isopropanol (Fig.14)
- Picral to develop difference between martensite and austenite:  
100 ml ethanol  
4 g picric acid  
10-60 sec. depending on carbon content.  
Rinse with water followed by isopropanol (Fig.13)

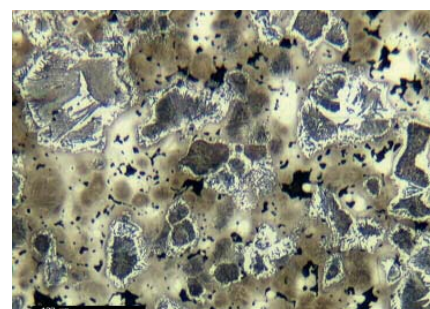


Fig.13: P/M steel with 0.5%C, diffusion alloyed with Ni, Cu and Mo. Etched in picral, showing areas of fine pearlite surrounded by ferrite, martensite, bainite and Ni-rich austenite.

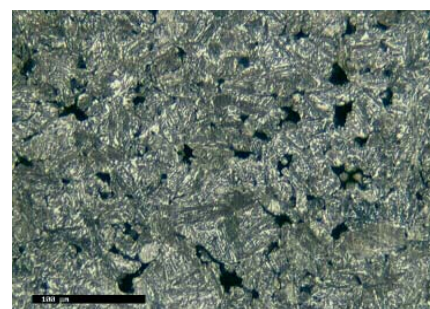


Fig.14: P/M steel with 0.8% C, pre alloyed with 1.5% Mo. Etched with nital, showing dense bainite



**Struers A/S**  
 Pederstrupvej 84  
 DK-2750 Ballerup, Denmark  
 Phone +45 44 600 800  
 Fax +45 44 600 801  
 struers@struers.dk

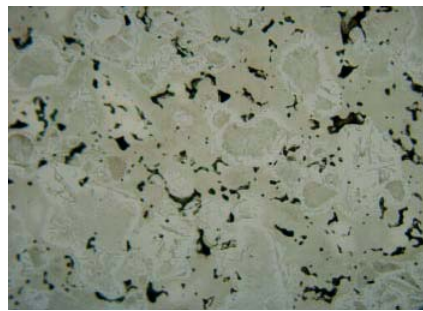


Fig.15: Etched too short

- For contrast between pearlite, bainite and martensite:  
 200 ml ethanol  
 4 g picric acid  
 1-2 ml nitric acid  
 20-100 sec, depending on carbon content and alloying elements.  
 Rinse with water followed by isopropanol

#### Stainless steels:

- Vilella's reagent:  
 45 ml glycerol  
 15 ml nitric acid  
 30 ml hydrochloric acid  
 30 sec to 5 min  
 Rinse thoroughly with water followed by isopropanol

#### Sintered tungsten carbides:

- Murakami's reagent  
 100 ml water  
 10 g potassium or sodium hydroxide  
 10 g potassium ferricyanide  
 Etching by immersion or swab etch  
 Rinse thoroughly with water followed by isopropanol

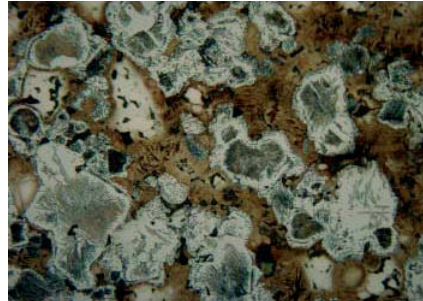


Fig.16: Etched too long, over etched

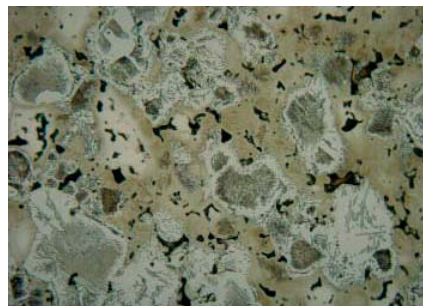


Fig.17: Correct etching

#### Summary

Powder metallurgy is a method of producing parts from metals which normally would not alloy easily, or from mixtures of metals and non-metals, taking advantage of their combined properties. Although the powder metallurgical process is costly, it has the advantage of economically producing large amounts of small and intricately shaped parts with homogenous structures. The density of compacted and sintered parts affects their strength, ductility and hardness. Therefore metallographic control of porosity is an integral part of quality control. During metallographic grinding, metal is pushed into the pores and if the subsequent polishing steps are not carried out properly, residual metal "lids" covering the pores will obstruct the evaluation of the correct porosity. Thorough grinding and polishing with diamond, and microscopic checks between polishing steps, assure a true representation of the structure. The given procedures for automatic preparation and the chemical etchants have been used successfully in day to day practical laboratory application and give reproducible results.

#### Author

Elisabeth Weidmann, Birgitte Nielsen,  
 Struers A/S, Copenhagen, Denmark  
 Judy Arner, Struers, Inc., Westlake, OH, USA

#### Acknowledgements

We wish to thank Höganäs AB, SE 26383 Höganäs, Sweden, for the co-operation and supplying information and material for reproduction: all photos of P/M components, diagram and SEM photos on page 2, Figs. 1, 2, 6, 7, 8, 9, 10, and 13-17. Special thanks go to M. Carlsson who generously shared her professional experience with us.  
 Figs. 4 and 5 courtesy Böhler Edelstahl GmbH, Kapfenberg, Austria

#### Bibliography

Höganäs Handbook for sintered components, 6. Metallography, Höganäs AB, 1999  
 Höganäs iron and steel powders for sintered components, Höganäs AB, 1998  
 Structure 24, Metallographic examination of P/M stainless steel matrix composites, N. Frees, Institute for Product Development, Technical University, Lyngby, Denmark, 1991  
 Metallographic etching, G. Petzow, ASM, Metals Park, Ohio, 1978



P/M steel with copper infiltration

**USA**  
**Struers Inc.**  
 24766 Detroit Road  
 Westlake, OH 44145-1598  
 Phone +1 440 871 0071  
 Fax +1 440 871 8188  
 info@struers.com

**CANADA**  
**Struers Ltd.**  
 7275 West Credit Avenue  
 Mississauga, Ontario L5N 5M9  
 Phone +1 905-814-8855  
 Fax +1 905-814-1440  
 info@struers.com

**SWEDEN**  
**Struers A/S**  
 Smältvägen 1  
 P.O. Box 11085  
 SE-161 11 Bromma  
 Telefon +46 (0)8 447 53 90  
 Telefax +46 (0)8 447 53 99  
 info@struers.dk

**FRANCE**  
**Struers S.A.S.**  
 370, rue du Marché Rollay  
 F- 94507 Champigny  
 sur Marne Cedex  
 Téléphone +33 1 5509 1430  
 Télécopie +33 1 5509 1449  
 struers@struers.fr

**NEDERLAND/BELGIE**  
**Struers GmbH Nederland**  
 Electraweg 5  
 NL-3144 CB Maassluis  
 Tel. +31 (0) 10 599 72 09  
 Fax +31 (0) 10 599 72 01  
 glen.van.vugt@struers.de

**BELGIQUE (Wallonie)**  
**Struers S. A. S.**  
 370, rue du Marché Rollay  
 F- 94507 Champigny  
 sur Marne Cedex  
 Téléphone +33 1 5509 1430  
 Télécopie +33 1 5509 1449  
 struers@struers.fr

**UNITED KINGDOM**  
**Struers Ltd.**  
 Unit 25a  
 Monkspath Business Park  
 Solihull  
 B90 4NZ  
 Phone +44 0121 745 8200  
 Fax +44 0121 733 6450  
 info@struers.co.uk

**IRELAND**  
**Struers Ltd.**  
 Unit 25a  
 Monkspath Business Park  
 Solihull  
 B90 4NZ  
 Phone +44 (0)121 745 8200  
 Fax +44 (0)121 733 6450  
 info@struers.co.uk

**JAPAN**  
**Marumoto Struers K.K.**  
 Takara 3rd Building  
 18-6, Higashi Ueno 1-chome  
 Taito-ku, Tokyo 110-0015  
 Phone +81 3 5688 2914  
 Fax +81 3 5688 2927  
 struers@struers.co.jp

**CHINA**  
 Struers Ltd.  
 Office 702 Hi-Shanghai  
 No. 970 Dalian Road  
 Shanghai 200092, P.R. China  
 Phone +86 (21) 5228 8811  
 Fax +86 (21) 5228 8821  
 struers.cn@struers.dk

**DEUTSCHLAND**  
**Struers GmbH**  
 Karl-Arnold-Strasse 13 B  
 D-47877 Willich  
 Telefon +49(02154) 486-0  
 Telefax +49(02154) 486-222  
 verkauf.struers@struers.de

**ÖSTERREICH**  
**Struers GmbH**  
 Zweigniederlassung Österreich  
 Ginzkeyplatz 10  
 A-5020 Salzburg  
 Telefon +43 662 625 711  
 Telefax +43 662 625 711 78  
 stefan.lintschinger@struers.de

**SCHWEIZ**  
**Struers GmbH**  
 Zweigniederlassung Schweiz  
 Weissenbrunnstrasse 41  
 CH-8903 Birmenndorf  
 Telefon +41 44 777 63 07  
 Telefax +41 44 777 63 09  
 rudolf.weber@struers.de

**CZECH REPUBLIC**  
**Struers GmbH**  
 Organizační složka  
 Havlíčkova 361  
 CZ-252 63 Roztoky u Prahy  
 Tel: +420 233 312 625  
 Fax: +420 233 312 640  
 david.cernicky@struers.de

**POLAND**  
**Struers Sp. z o.o.**  
 Oddział w Polsce  
 ul. Lirowa 27  
 PL-02-387 Warszawa  
 Tel. +48 22 824 52 80  
 Fax +48 22 882 06 43  
 grzegorz.uszynski@struers.de

**HUNGARY**  
**Struers GmbH**  
 Magyarországi fióktelep  
 Puskás Tivadar u. 4  
 H-2040 Budaörs  
 Phone +36 (23) 428-742  
 Fax +36 (23) 428-741  
 zoltan.kiss@struers.de

**SINGAPORE**  
**Struers A/S**  
 627A Aljunied Road,  
 #07-08 BizTech Centre  
 Singapore 389842  
 Phone +65 6299 2268  
 Fax +65 6299 2661  
 struers.sg@struers.dk