

ABPM 7100 Ambulatory Blood Pressure Monitor



Directions for Use

WelchAllyn[®]

Advancing Frontline Care[™]

© 2014 Welch Allyn, Inc. To support the intended use of the product described in this publication, the purchaser of the product is permitted to copy this publication for internal distribution only, from the media provided by Welch Allyn.

Caution: Federal US law restricts sale of the device identified in this manual to, or on the order of, licensed physicians.

Welch Allyn does not accept liability for injuries or unlawful or improper use of the product which may result from the fact that the product is not used in accordance with the instructions, cautions and warnings, as well as the indications for use published in this manual.

Welch Allyn is a registered trademark of Welch Allyn, Inc., ABPM 7100 and CardioPerfect are trademarks of Welch Allyn, Inc.

The copyright for the firmware in this product remains with the manufacturer of this device. All rights reserved. The firmware may not be read out, copied, decompiled, redeveloped, disassembled or brought into any human-readable format. This does not pertain to the sales of firmware or a firmware copy. All usage and ownership rights to the Software remain with I.E.M. GmbH.

Welch Allyn Technical Support:

<http://www.welchallyn.com/about/company/locations.htm>



I.E.M. GmbH
Cockerillstr. 69
52222 Stolberg
Germany

Manufactured for Welch Allyn










Table of contents

Symbols	4
Introduction	5
Preliminary note	5
About these directions for use	5
Clinical trials	5
CE Mark	5
Content	6
Direction for use	6
Intended use	6
Indications for use	6
Improper use	6
Essential Performance	7
Side effects of 24-hour blood pressure monitoring	7
Product description	7
Introduction	7
The ABPM 7100	8
Technical Data	10
Preparing the ABPM 7100	11
Safety instructions	11
Inserting the batteries	12
Activating the device	13
Setting the time / date	13
Clearing the memory	13
Transferring patient data (ID)	14
Setting measurement logs	14
Selecting a suitable cuff	15
Applying the ABP Monitor and cuff	15
Connecting the cuff tubing to the ABPM 7100	17
Positioning the patient for measurement	18
Measurement process	18
Safety instructions	18
Initial measurement	20
24-hour measurement	20
Performing a measurement	20
Cancelling a measurement	21
Unsuccessful measurement	21
Care and Maintenance	21
Cleaning and disinfection	21
Maintenance plan	22
Troubleshooting	23
Basic error sources	23
Transmission error	23
Checklist	23
Error codes	24
Limited Warranty	27
Service Policy	27
EMC Guidelines and Manufacturer Declaration	28
Patient Information - operation of the ABPM 7100	31


4 - Symbols

Symbols


Documentation symbols

	WARNING The warning statement identifies an immediate threat. Non-adherence may lead to the most severe injuries and to death		CAUTION The caution statement identifies a possible hazard. Non-adherence may lead to minor or moderate injuries
Attention	The attention statement marks possible material damage. Non-adherence may lead to damage to the device or its accessories	Note	The note statement marks further information on the ABPM 7100 or its accessories
	INTERNAL REFERENCE Marks references within the document to further information		EXTERNAL REFERENCE Marks references to external documents containing further optional information
	Mandatory – Consult Directions for Use		Meets essential requirements of European Medical Device Directive 93/42/EEC
	Consult Directions for Use, Electronic version available at Welchallyn.com , or Hard copy DFU available from Welch Allyn within 7 days.		



Power symbols

	Battery symbol indicates the type of power supply		Nonionizing electromagnetic radiation
---	---	---	---------------------------------------










Connectivity symbols

	FCC ID and IC		Bluetooth Connectivity
---	---------------	---	------------------------

Shipping, storing and environment symbols

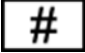	Separate the device from other disposables for recycling.
	See www.welchallyn.com/weee

Miscellaneous symbols

	Manufacturer		Date of manufacture yyyy-mm
	Reference/Model number		Serial number
	Reorder/Catalog number		Batch code
	Global Trade Item Number		Protection class
	Defibrillation-proof type BF applied part		

Introduction

Preliminary note

With the ABPM 7100 24-hour blood pressure measuring device, you now have an Ambulatory Blood Pressure Monitoring System (ABPM System) at your disposal.

The ABPM 7100, also specified as ABP Monitor, can be prepared for a new patient in just a few minutes. This permits the optimum use of the ABP Monitor and allows you to process one 24-hour profile per day.

The ABPM 7100 can therefore be quickly integrated into everyday practice life. The recorded blood pressure values must be evaluated with the intended software.

In combination with the Hypertension Management Software and an appropriate license, the ABPM 7100 is also able to process a haemodynamic analysis of the recorded pulse waves.

About these directions for use

These directions for use will familiarize you with the use of the ABPM 7100 and its accessories.

The software **CardioPerfect Workstation (CPWS)** is used for measurement evaluation.

Upgrades for haemodynamic evaluation may also be purchased from Welch Allyn. Please contact Welch Allyn for further information.

With reference to specific version characteristics, only the parts relevant for your respective version will apply.



- Please refer to the CPWS directions for use for software operating instructions.
- For the upgrades please refer to the respective directions for use to operate the Hypertension Management Software (HMS), version 5.0 and above.

Note

These directions for use explain the ABPM 7100 and its accessories in the sequence in which you setup the device for a blood pressure measurement, followed by the installation, initial operation, measurement preparation, placement on the patient and the evaluation. Individual functions are only explained when they are needed. You will therefore be familiarized with the ABPM 7100 on a step-by-step basis.

This directions for use must be kept with the product for later use!

Clinical trials

The blood pressure measuring device ABPM 7100 fulfills the requirements of the ESH (European Society of Hypertension).

The blood pressure measuring device ABPM 7100 fulfills the requirements of the BHS (British Hypertension Society).

CE Mark



The ABPM 7100 fulfills the requirements of the following directives:

- Directive 93/42/EEC (MDD)
- Directive 1999/5/EC (R&TTE)
- Directive 2011/65/EU (RoHS)

The ABPM 7100 bears the CE mark.

6 - Direction for use

Content

Standard

1. ABPM 7100 Monitor
2. Pressure Cuff – Size “Adult”
3. Carrying Pouch
4. PC Interface Cable
5. 4x AA Alkaline Batteries
6. ABPM 7100 Directions For Use
7. Calibration Notice
8. CPWS Software (*dependent on set*)
9. Pressure Cuff – Size “Adult Plus” (*dependent on set*)

HMS Option

1. HMS Software
2. Bluetooth®-Dongle
3. Quick Start Guide (*dependent on upgrade option*)
4. Version dependent 16 digit License Code (*dependent on upgrade option*)

Warning

Risk of injury posed by the use of other accessories. The use of unapproved accessories may lead to incorrect measurement results.

- Only use accessories approved and distributed by the manufacturer and Welch Allyn.
- Check the accessories regarding the manufacturer’s information before first use.

Direction for use

Intended use

The ABPM 7100 is intended for clarifying the blood pressure status and for use as a diagnostic aid for an individual adult patient (in the patient’s environment). The ABPM 7100 may only be used under medical supervision and after detailed instruction has been provided by the doctors or healthcare professionals. The ABPM 7100 in combination with the Hypertension Management Software (HMS) provides a derived ascending aortic blood pressure wave form and a range of central arterial indices.

Indications for use

- The ABPM 7100 is an automated microprocessor controlled ambulatory blood pressure monitor (ABPM) which monitors, accumulates and stores: heart beat (rate), systolic and diastolic data of an individual adult patient (in the patient’s environment) for a session which may last 24 hours.
- The ABPM 7100 in combination with the Hypertension Management Software (HMS) provides a derived ascending aortic blood pressure wave form and a range of central arterial indices.
- It is used with a standard cuff for blood pressure measurement.
- It is used in those patients where information related to the ascending aortic blood pressure is desired but in the opinion of the physician, the risk of cardiac catheterization procedure or other invasive monitoring may outweigh the benefits.

Improper use

- The ABPM 7100 may not be used for any purpose other than for the blood pressure monitoring procedures described here!
- Due to the strangulation risk posed by tubing and cuff, the ABPM 7100 may not be used on patients with limited cognitive abilities, patients under anesthetics and must not be placed within reach of unsupervised children.
- The ABPM 7100 must not be used on neonates and children under the age of 3 years. The device has not been tested with pregnant women, including preeclamptic patients.
- The ABPM 7100 must not be used for blood pressure monitoring purposes during surgeries.
- The ABPM 7100 is not intended for alarm-triggering monitoring purposes in intensive care units.

Essential Performance

The main performance features are defined as blood pressure measurement with:

- Error tolerances of the pressure gauge and measurement results are within required limits (IEC 80601-2-30).
- Maximum change value in blood pressure determination is as specified in IEC 80601-2-30.
- Cuff pressurization remains within specified limits (IEC 80601-2-30).
- An error is issued in the event that successful blood pressure measurement is impossible.

The ABPM 7100 does not issue ALARMS pursuant to IEC 60601-1-8 and is not intended for use in connection with HF surgical equipment or to clinically monitor patients in intensive care units.

Basic safety means that the patient cannot be endangered by any automatic device procedure. During any unclear conditions, the ABPM 7100 must therefore transfer into the safe **Standby** mode, during which the ABPM 7100 cannot automatically inflate the cuff, while this can be manually triggered by pushing the **START** button.

In this context, any interruption of a measurement or in automatic operation by an external influence, or the ability of the ABPM 7100 to test error conditions, is regarded as the retention or restoration of basic safety, and not as non-adherence to the main performance features.

Side effects of 24-hour blood pressure monitoring

- 24-hour blood pressure monitoring is a commonly practiced, accepted and valued measurement technique, used in daily diagnostics and treatment monitoring.
- When considering 24-hour blood pressure monitoring, always check whether the patient displays coagulation disorders or is being treated with anticoagulants. As with occasional blood pressure measurements, petechial haemorrhage (bleeding) may occur on the measured arm despite a correctly seated cuff. The innate patient-dependent risk resulting from treatment with anticoagulants or in patients with coagulations disorders arises irrespective of the type of monitoring device.

Product description

Introduction

The ABPM 7100 System consists of two main components:

- The ABPM 7100 with cuffs and accessories
- Patient management software for the doctor to evaluate the measurement results

With the **CPWS** software the ABPM 7100 can be prepared for measurement, transfer stored measurement results to the PC, display transferred measurements on the screen in various formats such as graphics, lists and statistics and print out measurement results. Optional is the possibility to evaluate the measurement results with the **HMS** and its upgrades.

The ABPM 7100 can be prepared immediately for the next patient. With little practice this procedure can be completed in just a few minutes. This allows the doctor to use the ABPM 7100 around the clock on each work day.

The ABPM 7100 is designed to allow recording and display of a blood pressure profile throughout the day and at night. Additional parameters such as nocturnal values and blood pressure fluctuations are recognized. This permits the doctor to prescribe optimal medical treatment for each individual.

Measurement with the ABPM 7100 can be either automated or be manually controlled by the user. In order to start a series of automatic measurement, the user must initiate the first measurement by pressing the **START** button and the doctor should check the reliability of the first measurement.

During the first measurement, the cuff is inflated in increments, to determine the cuff pressure required to measure the systolic blood pressure value. The maximum required inflation pressure is stored and applied by direct inflation during the subsequent automatic measurements. This procedure is called **AFL – Auto Feedback Logic**.

8 - Product description

The ABPM 7100

Components

- 1 Cuff connection
- 2 ON/OFF button
- 3 LCD-Display
- 4 START button
- 5 DAY/NIGHT button
- 6 EVENT button
- 7 PC Interface cable port



The Buttons



ON/OFF

The **ON/OFF** button turns the ABPM 7100 on and off. To prevent unintended activation, the ABPM 7100 turns on or turns off only when the button is pressed for more than 2 seconds.



START

The **START** button serves to

- initiate a manual measurement to ascertain whether the ABPM 7100 is working correctly.
- initiate a 24-hour measurement.
- perform a measurement outside the specified measurement cycle.



DAY/NIGHT

The **DAY/NIGHT** button is used to differentiate between waking and sleeping phases during the measurement, which is important for statistics and the graphic displays.

The patient is instructed to press the **DAY/NIGHT** button upon going to bed and again, when getting up in the morning. This individually adapts the measurement interval to the patient and assists you in the analysis of the blood pressure profile.



EVENT

The patient uses the **EVENT** button to document the time of medication or to record any events which may cause the blood pressure to rise or fall. Pressing the button will trigger a measurement, the patient should note the reason for pressing the **EVENT** button in the event log.

LCD Display

The LCD display is located on the front of the ABPM 7100 casing. It displays useful information for the doctor and the patient regarding measurement data, monitor settings and measurement errors. When the **START** button is pressed, the number of previously registered measurements will be shown before starting a manual measurement.

Audible signals

Individual or multiple beeps of audible signals are used. The following table explains the meaning of the beeps:

1 beep	<ul style="list-style-type: none"> • Switching ON/OFF • Starting and ending a measurement (except at night intervals) • Removal of the interface cable • Establishing and ending Bluetooth® communication • Measurement errors
3 beeps	<ul style="list-style-type: none"> • System errors
Continuous beeps	<ul style="list-style-type: none"> • Severe system errors (e.g. cuff pressure is higher than 15 mmHg for longer than 10 seconds outside the measurement)
Combined beeps	<ul style="list-style-type: none"> • Manual deletion of measurement, 1 beep followed by 5 beeps 2 seconds later

Cuff connection

- The cuff connection is located at the top of the ABPM 7100 casing.
- The cuff is connected to the ABPM 7100 via a metal connector.


Attention

Measurement errors


- The cuff connection must always engage with an audible “CLICK”. A poor connection between the ABPM 7100 and cuff will result in measurement errors.

The Arm Cuff

- 1 The arm cuff
- 2 Air tube
- 3 Air tube connection

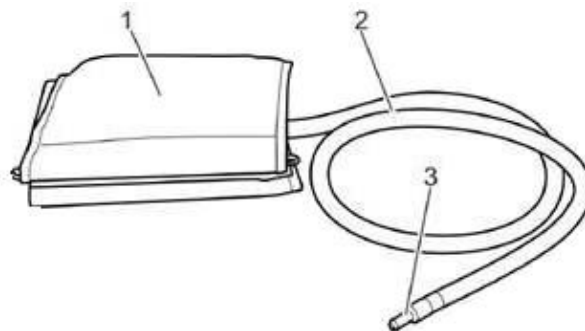


PC Interface Cable

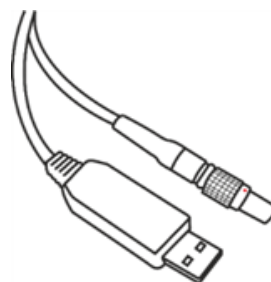

In order to read data from the ABPM 7100, the interface cable must be connected to an USB slot on a PC.

10 - Product description

PC Interface Cable Port

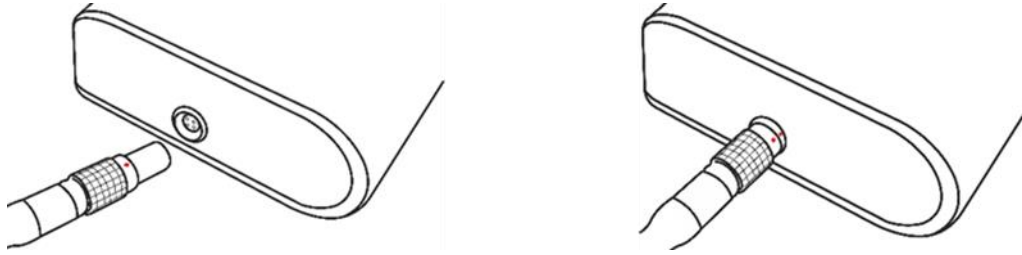

- The connecting port for the PC interface cable is located at the bottom of the ABPM 7100 casing.
- The red dot on the plug must align with the red dot on the port before plugging.
- To disconnect, pull on the knurled ring of the connector.

Connecting the ABPM 7100 to the PC

To transfer the data from the ABPM 7100, ensure that the interface cable is connected correctly to an USB port on the PC and to the interface cable port on the device.

Technical Data

Measurement pressure range:	Systolic 60 to 290 mmHg Diastolic 30 to 195 mmHg
Accuracy:	+/- 3 mmHg in display range
Static pressure range:	0 to 300 mmHg
Pulse range:	30 to 240 beats per minute
Procedure:	oscillometric
Measurement intervals:	0, 1, 2, 3, 4, 5, 6, 10, 12, 15, 20 or 30 measurements per hour
Measurement logs:	4 adjustable interval groups
Memory capacity:	300 Measurements
Battery capacity:	> 300 Measurements
Operating temperatures:	+5°C to +40°C
Operating humidity:	15% to 93%
Storage environment:	-25°C to +70°C and 15% to 93% humidity
Dimensions:	121 x 80 x 33 mm
Weight:	approx. 220 g excluding batteries
Power supply:	2 Ni-MH batteries with 1,2 V each and min. 1500 mAh (AA, Mignon) or 2 Alkali 1,5 V batteries (AA, Mignon, LR6)
Interfaces:	USB-Interface cable Bluetooth® (Class 1 / 100 m and max. 100 mW with 2,402 GHz to 2,480 GHz) only available with optional HMS software
Expected operational device life:	5 Years
Expected operational cuff life:	6 Months

Preparing the ABPM 7100

Safety instructions

Warning

Risk of strangulation posed by the shoulder strap and cuff tubing.

- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Do not place the shoulder strap and cuff tubing around the patient's neck.
- Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.
- Measurement can be interrupted at any stage by pushing any of the buttons. This deflates the cuff and the device can be removed.

Caution

Risk of injury caused by incorrect application of the device.

- The doctor must ensure that, due to the patient's medical condition, the use of the device does not result in impaired blood circulation in the arm.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- While it is still attached to a patient, the device may never be connected to a PC or other device.
- Instruct the patient to place the device in such a way that, while the cuff is inflated, the tubing is not compressed or kinked, especially during sleep.
- Petechiae, haemorrhages or subcutaneous haematoma may occur in some patients.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Caution

Risk of injury caused by incorrect application of the cuff.

- The doctor must ensure that, due to the patient's medical condition, the use of the device does not result in impaired blood circulation in the arm.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- It is imperative that you instruct the patient in the correct seating of the cuff.
- Inform the patient that the cuff may only be used on the upper arm.
- Ensure that neither the shoulder strap nor the cuff tubing can ever wrap around the patient's neck. Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to place the device in such a way that, while the cuff is inflated, the tubing is not compressed or kinked, especially during sleep.
- Petechiae, haemorrhages or subcutaneous haematoma may occur in some patients.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Caution

Intolerances caused by the use of disinfectants.

- Wash to remove residues.
- Wash the cuff sleeve with a mild detergent in the washing machine at max. 30°C without spinning.

12 - Preparing the ABPM 7100

Inserting the batteries

Attention

Damage to device

- Do not wear the ABPM 7100 while showering! If you suspect that liquid has entered the device while cleaning or using it, the device shall no longer be used on the patient.
- If the device was exposed to moisture, switch off the device and remove the batteries.
- Inform your service immediately and send the device in for inspection.
- The device may not be operated around MRI scanners or in the immediate vicinity of other medical electrical equipment.
- During a defibrillator discharge, the device shall not be in contact with the patient. Such a discharge may damage the ABPM 7100 and cause it to display incorrect values.
- The ABPM 7100 is not suitable for simultaneous use with HF surgical equipment.
- Measurement can be interrupted at any stage by pushing a random button. This deflates the cuff and the device can be removed.

Attention


Device function

- Although zinc-carbon batteries may indicate sufficient voltage during a battery test, their output is frequently insufficient to perform 24-hour measurements.
- To conduct measurements over 48 hours, you will need 2 additional batteries as replacement after 24 hours.

Attention

Damage to device

- Do not open the casing. Once the device is opened, all warranties will lapse.

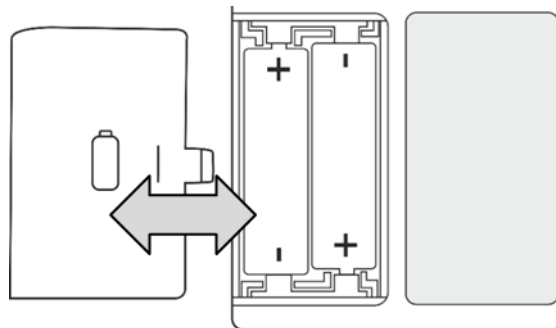

Open the batteries compartment on the rear of the ABPM 7100 casing to insert the batteries into the ABPM 7100 according to the battery polarities (+ / -) and close the compartment.

Note

- Always use fully charged batteries for a new measurement
- Only use undamaged batteries
- Please remove the batteries if the device has not been used for a longer period
- When inserting the batteries, please ensure correct polarity

Activating the device

Attention

Damage to device

- Do not wear the ABPM 7100 while showering. If you suspect that liquid has entered the device while cleaning or using it, the device shall no longer be used on the patient.
- In the device was exposed to moisture, switch off the device and remove the batteries.
- Inform your service immediately and send the device in for inspection.
- The device may not be operated around MRI scanners or in the immediate vicinity of other medical electrical equipment.
- During a defibrillator discharge, the device shall not be in contact with the patient. Such a discharge may damage the ABPM 7100 and cause it to display incorrect values.
- The ABPM 7100 is not suitable for simultaneous use with HF surgical equipment.
- Measurement can be interrupted at any stage by pushing a random button. This deflates the cuff and the device can be removed.

Attention

Hygiene

- Ensure hygiene in accordance with the maintenance schedule.

Always check the condition of the ABPM 7100 by observing the initial display shown on the device shortly after turning it on and before handing it to a patient. The ABPM 7100 performs a self-test. In addition, a beep sounds to check the speaker. The following should be displayed in this sequence:

Test	Display	Comment
Battery condition (volts)	2.85	(At least 2.75 volts for NiMH batteries and at least 3.10 volts for alkaline batteries)
Display Segment Test	999:999 to 000:000	The display of the figures (999:999 to 000:000) is accompanied by all other symbols of the LCD in succession. Check whether all segments are correctly and fully displayed (the complete program code is checked for correctness in the background)
Current 24h time	21:45	hh:mm

If the internal test detects an error, the ABPM 7100 will indicate “**E004**” on the display and emit an audible signal. For safety reasons, the use of the ABPM 7100 will be locked. The faulty ABPM 7100 unit should be send back immediately for repairs to your dealer or to Welch Allyn.

Setting the time / date

The ABPM 7100 has an internal buffer battery allowing the time to continue even if the batteries have been removed. Nevertheless the time and the date should be checked before every measurement series.

The time and date can be set automatically with the patient management software.

Alternatively the time and date can be set manually. Press and hold the **START** button and then press the **EVENT** buttons to enter the **Set Time** mode. Use the **START** button to select the appropriate item and use the **EVENT** button to jump to the next display item.

Clearing the memory

The device memory must be cleared before every measurement series, i.e. blood pressure data from the previous patient must not remain in the memory.

If there are existing data, the memory can be cleared with the delete function of the analysis software.

Alternatively the data can be cleared manually. Press and hold the **START** button for a minimum of 5 seconds until “**cLr**” is displayed. Within the next 5 seconds press and hold the **EVENT** button for at least 2 seconds to confirm the deletion of the stored measurements. The device emits a single beep to indicate that the memory is cleared.

14 - Preparing the ABPM 7100

Transferring patient data (ID)

The ABPM 7100 must be prepared by transferring patient data (ID) with the help of the patient management software, so that correct data allocation is possible when it is read out after measurement.

Please refer to the respective patient management software manual for how to transfer patient data (ID) to the ABPM 7100.

Setting measurement logs

In the patient management software you can optionally choose between eleven (1-11) logs. A log serves to set the measurement intervals. As soon as you have conducted a measurement, the log can only be changed once you have fully deleted all data.

Manual log settings

For manual log setting, press and hold the **DAY/NIGHT** button while simultaneously pressing the **EVENT** button. Use the **START** button to change the log and confirm with the **EVENT** button.

Log	Day-Time	Night-Time	Measurements per hour	Audible signal	Display of measured values
1	08:00	23:59	4	YES	YES
	00:00	07:59	2	NO	
2	08:00	22:59	4	YES	YES
	23:00	07:59	1	NO	
3	07:00	21:59	4	YES	NO
	22:00	06:59	2	NO	
4	08:00	23:59	4	YES	NO
	00:00	07:59	2	NO	
5	18:00	09:59	4	YES	YES
	10:00	17:59	2	NO	
6	07:00	23:59	4	YES	YES
	00:00	06:59	2	NO	
7	06:00	22:59	4	YES	NO
	23:00	05:59	2	NO	
8	07:00	08:59	6	YES	YES
	09:00	23:59	4	YES	
	00:00	06:59	2	NO	
9	09:00	08:59	30	NO	YES
10	08:00	07:59	30	YES	NO
11	08:00	23:59	4	YES	YES
	00:00	07:59	2	NO	

Setting the logs via software

To set the logs via software please refer to the respective patient management software manual.

- Note**
- Logs 1, 2, 10 and 11 are set by default but can be changed via the patient management software.
 - Log 5 is suitable for nighttime activities (night shift).
 - Log 9 is designated as "Schellong-Test".
 - Log 10 automatically sends the measurement values to your doctor's PC via Bluetooth®. Bluetooth® communication is not supported with the CPWS software.
 - Log 11 is only available to upgraded ABPM 7100 systems in connection with HMS from version 5.0. Blood pressure measurement intervals and the 24h PWA can be set separately here. Please contact Welch Allyn for further information.

Selecting a suitable cuff

Caution

Risk of injury caused by incorrect application of the cuff.

- The doctor must ensure that, due to the patient's medical condition, the use of the device does not result in impaired blood circulation in the arm.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- It is imperative that you instruct the patient in the correct seating of the cuff.
- Inform the patient that the cuff may only be used on the upper arm.
- Ensure that neither the shoulder strap nor the cuff tubing can ever wrap around the patient's neck. Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to place the device in such a way that, while the cuff is inflated, the tubing is not compressed or kinked, especially during sleep.
- Petechiae, haemorrhages or subcutaneous haematoma may occur in some patients.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Caution

Intolerances caused by the use of disinfectants.

- Wash to remove residues.
- Wash the cuff sleeve with a mild detergent in the washing machine at max. 30°C without spinning.

The correct cuff size is important for correct blood pressure measurement. To obtain reproducible measurements, standardized measuring conditions are needed. Measure the circumference of the upper arm and select the appropriate cuff:

Welch Allyn Size Number	Upper Arm Circumference	Cuff
10	20 – 24 cm	Small Adult
11	24 – 32 cm	Adult
11L	32 – 38 cm	Adult Plus
12	38 – 55 cm	Large Adult

Applying the ABP Monitor and cuff

Warning

Risk of strangulation posed by the shoulder strap and cuff tubing.

- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Do not place the shoulder strap and cuff tubing around the patient's neck.
- Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.
- Measurement can be interrupted at any stage by pushing any of the buttons. This automatically deflates the cuff and the device can be removed.

Warning

Poor circulation caused by continuous cuff pressure.

- Do not kink the connecting tubing.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Ensure the correct placement of the shoulder strap and cuff tubing.
- Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

16 - Preparing the ABPM 7100

Warning

Placement and inflation of the cuff over a wound may lead to further injuries.

Placement and inflation of the cuff on any limb with an intravascular access or under intravascular treatment or an arteriovenous (A-V) shunt may lead to temporary interruption of circulation and therefore to further patient injury.

Placement and inflation of the cuff on the arm at the side of a breast amputation may lead to further injury.

- Examine the patient for wounds, bandages, etc.
- Question the patient regarding previous treatments.
- Observe the patient closely.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Caution

Intolerances caused by the use of disinfectants.

- Wash to remove residues.
- Wash the cuff sleeve with a mild detergent in the washing machine at max. 30°C without spinning.

Caution


Risk of injury caused by incorrect application of the cuff.

- The doctor must ensure that, due to the patient's medical condition, the use of the device does not result in impaired blood circulation in the arm.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- It is imperative that you instruct the patient in the correct seating of the cuff.
- Inform the patient that the cuff may only be used on the upper arm.
- Ensure that neither the shoulder strap nor the cuff tubing can ever wrap around the patient's neck. Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to place the device in such a way that, while the cuff is inflated, the tubing is not compressed or kinked, especially during sleep.
- Petechiae, haemorrhages or subcutaneous haematoma may occur in some patients.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Caution

Risk of injury caused by incorrect application of the device.

- The doctor must ensure that, due to the patient's medical condition, the use of the device does not result in impaired blood circulation in the arm.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- While it is still attached to a patient, the device may never be connected to a PC or other device.
- Instruct the patient to place the device in such a way that, while the cuff is inflated, the tubing is not compressed or kinked, especially during sleep.
- Petechiae, haemorrhages or subcutaneous haematoma may occur in some patients.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

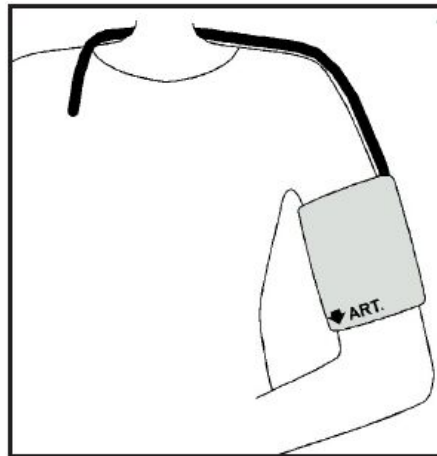

Applying the ABP Monitor and cuff:

1. Position the carrying pouch on the right side of the patient. By varying the length of the pouch strap, it can be worn around the hips or around the shoulders.
2. Alternatively a normal belt matching the clothes can be used.
3. Fit the cuff onto the patient.
The correct cuff seating is very important for correct blood pressure measurement.
4. Align the cuff so that no part of the cuff tubing is kinked. In this regard, the tube connection on the cuff must face upwards.
5. Align the cuff so that the lower edge of the cuff is approximately 2 cm above the inside of the patient's elbow.
6. Tighten the cuff around the upper arm until one finger can be introduced under the cuff.
7. It is imperative that the artery symbol is positioned on the brachial artery. If you have aligned the cuff correctly, the metal bar will lie on the outside of the upper arm (on the elbow side), whereby the cuff sleeve must cover the skin under the metal bar.
8. Guide the tubing through the shirt's row of buttons and out of the clothing, behind the nape of the neck to the ABPM 7100 on the right side of the body.
9. The cuff can be worn on the naked upper arm or over a thin shirt sleeve.
10. The placement of the pressure tube must guarantee the upper arm's free motion.

Connecting the cuff tubing to the ABPM 7100

1. Push the tube firmly onto the connection, with the cuff tubing engaging with an audible "CLICK" (to detach, simply pull back the knurled ring).
2. Before measurement, check to ensure that tubing, ABPM 7100 and cuff are positioned correctly. The ABPM 7100 is ready for measurement only once this is ensured.

18 - Measurement process

Positioning the patient for measurement

The patient should take the following position during blood pressure measurement:

- Comfortably seated
- Legs uncrossed
- Feet flat on the ground
- With support to the back and the arms
- With the cuff center at one level with the right atrium

- Note**
- During measurement, the patient should be as relaxed as possible and may not speak unless he wants to report any discomfort!
 - Allow for 5 minutes of relaxation before recording the first measurement value.
 - Blood pressure measurements can be influenced by the patient's position (standing, sitting, lying), by exertion or the patient's physiological state. Exclude these influencing factors to the greatest possible degree!

Measurement process

Safety instructions

Warning

Risk of strangulation posed by the shoulder strap and cuff tubing.

- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Do not place the shoulder strap and cuff tubing around the patient's neck.
- Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.
- Measurement can be interrupted at any stage by pushing any of the buttons. This automatically deflates the cuff and the device can be removed.

Warning

Poor circulation caused by continuous cuff pressure.

- Do not kink the connecting tubing.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Ensure the correct placement of the shoulder strap and cuff tubing.
- Always place the cuff tubing under the outer clothing (even at night).
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Warning

Poor circulation due to overly frequent measurements.

- Check the date of the last measurement.
- Inform the patient about this warning.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Observe the patient closely.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

⚠ Warning

If the patient is wearing an additional ME device on the same limb for monitoring purposes, the placement and inflation of the cuff may trigger the temporary loss of the existing ME device's function.

The operation and use of the automated non-invasive blood pressure monitoring device may lead to longer impaired blood circulation in the patient or respective limb.

- Examine the patient.
- Question the patient regarding previous treatments.
- Observe the patient closely.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

⚠ Caution

Risk of injury caused by incorrect application of the cuff.

- The doctor must ensure that, due to the patient's medical condition, the use of the device does not result in impaired blood circulation in the arm.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- It is imperative that you instruct the patient in the correct seating of the cuff.
- Inform the patient that the cuff may only be used on the upper arm.
- Ensure that neither the shoulder strap nor the cuff tubing can ever wrap around the patient's neck. Always place the cuff tubing and the outer clothing (even at night).
- Instruct the patient to place the device in such a way that, while the cuff is inflated, the tubing is not compressed or kinked, especially during sleep.
- Petechiae, haemorrhages or subcutaneous haematoma may occur in some patients.
- Instruct the patient to switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

⚠ Caution

Intolerances caused by the use of disinfectants.

- Wash to remove residues.
- Wash the cuff sleeve with a mild detergent in the washing machine at max. 30°C without spinning.

Attention

Damage to device

- Do not wear the ABPM 7100 while showering. If you suspect that liquid has entered the device while cleaning or using it, the device shall no longer be used on the patient.
- In the device was exposed to moisture, switch off the device and remove the batteries.
- Inform your service immediately and send the device in for inspection.
- The device may not be operated around MRI scanners or in the immediate vicinity of other medical electrical equipment.
- During a defibrillator discharge, the device shall not be in contact with the patient. Such a discharge may damage the ABPM 7100 and cause it to display incorrect values.
- The ABPM 7100 is not suitable for simultaneous use with HF surgical equipment.
- Measurement can be interrupted at any stage by pushing a random button. This deflates the cuff and the device can be removed.

Attention

Hygiene

- Ensure hygiene in accordance with the maintenance schedule.

20 - Measurement process

Attention

Measurement errors

- The use of components not included in the scope of delivery may lead to measurement errors. Only use the accessories offered by Welch Allyn.
- Although the ABPM 7100 fulfills all EMC standard requirements, it should not be exposed to any strong electromagnetic fields, as this may lead to malfunctions outside the limit values.
- Medical electrical equipment is subject to special EMC precautions. Please observe the directives attached.
- The cuff tubing between the ABPM 7100 and the cuff may not be knotted, compressed or pulled apart.
- The cuff connection must always engage with an audible “CLICK”. A loose connection between the tubing and the device leads to measurement errors.

Note

- Severe malfunctions are indicated by a continuous beep.
- In the event of a continuous beep, switch off the device, remove the cuff and inform your doctor.
- Hand the data sheet “Patient information – operation of the ABPM 7100” to each patient. The data sheet is attached as a copy template.
- Portable and mobile HF communication equipment may influence medical electrical devices.
- Extreme temperatures, humidity or air pressure can influence measurement accuracy. Please observe the operating conditions.
- The pulse wave analysis provides additional indicators for possible risks, but is not permissible as a sufficient indicator for individual illness or as a treatment recommendation.

Initial measurement

Note

An initial measurement is required for starting the measurement log. The initial measurement must be checked by a physician for plausibility!

24-hour measurement

1. Ensure sufficient battery voltage. At least 2.6 volts for NiMH batteries and at least 3.10 volts for alkaline batteries!
2. The doctor must go through these instructions together with the patient before a 24-hour measurement.
3. The doctor must explain the possible hazards in detail on the basis of the warnings above!
4. Ensure that the patient has understood all functions and observable points!



Safety:

For your own safety during the following steps, please observe the safety instructions at the start of this chapter, as well as the functional overview.

Performing a measurement

1. To trigger a measurement, press the **START** button.
 - The number of stored measurements will be shown on the LCD display.
 - An audio beep will announce the upcoming measurement.
 - Manual measurement will start.
2. The patient should stay calm during the measurement process, until the measurement is completed. Allow your arm to hang loose, or place your lower arm loosely on the table or on a support whilst sitting. Avoid any movement!
3. Doctor: Please check the values of the first measurement for plausibility, so that subsequent automatic measurements can be processed correctly and correct cuff position is ensured.
4. In the event of an error measurement, please follow the instructions in sections **Measurement preparations** and **Troubleshooting**.

Cancelling a measurement

A measurement will be cancelled by pressing any buttons during the measurement process. The LCD display will then show **-StoP-** and the ABPM 7100 will beep 5 times. This cancellation will be stored in the measurement value table under **Cancel**.

Unsuccessful measurement

1. If the display shows errors, reexamine the correct procedure during set-up and positioning of the device.
2. Dismiss the patient only after a successful manual measurement!
Inform the patient sufficiently in order to explain the situation!
3. Repeat the measurement.
4. If the display still shows errors, repeat the initial operation process.
5. For further troubleshooting measures and faults removal, please refer to the **Troubleshooting** section.

- Note**
- Severe malfunctions are indicated by a continuous beep.
 - In the event of a continuous beep, switch off the device, remove the cuff and inform your doctor.

Care and Maintenance

To ensure the optimal functionality of the ABPM 7100 regular care and maintenance of the unit is required.

Attention

Damage to device

- Do not open the casing. Once the device is opened, all warranties will lapse.

Cleaning and disinfection

The user (doctor) decides whether and when the cuff sleeve should be disinfected for hygienic reasons (e.g. after every use).

The following agents are recommended for disinfecting the cuff sleeve:

- Terralin Liquid (Manufacturer: Schülke & Mayr)
- Promanum N (Manufacturer: B. Braun)

The use of disinfectants not recommended by Welch Allyn shall render the user responsible for proof of safe application.

- Note** It is imperative that you observe the manufacturer's information regarding the use of these products. Allow the agents to dry off completely.

Caution

Intolerances caused by the use of disinfectants: Some patients display intolerance (e.g. allergies) to disinfectants or their components.

- Never use disinfectants which leave residues on the product or which are not suitable for contact with the skin.
- Carefully wash the cuff to remove residues.

22 - Care and Maintenance

Cleaning the ABP Monitor and carrying pouch

Attention

Damage to the ABP Monitor and carrying pouch caused by the use of solvents

- Only use a moistened cotton cloth to clean the ABPM 7100 Monitor.
- Do not use strong or solvent-based additives.
- Ensure that no liquid enters the device.
- If liquid does penetrate the device, switch it off immediately and return it to your Welch Allyn specialist for inspection.

1. Use only a moistened cotton cloth to clean the ABPM 7100.
2. Read the safety instructions carefully and observe them closely during cleaning.

Cleaning the cuff sleeve, bladder and tubing

Attention

Damage to the cuff caused by disinfectants

- Do not submerge the cuff sleeve in disinfectants.
- Avoid disinfecting the cuff bladder and connected rubber tubing.

Attention

Damage to the cuff sleeve during washing

- Always close the Velcro strip before washing!
- Wash the cuff sleeve with a mild detergent in the washing machine at max. 30°C. Do not spin.
- Do not use fabric softeners or other washing aids (e.g. hygiene rinses, textile deodorants). These agents can leave behind residues and damage the material.
- The cuff sleeve is not suitable for drying in a dryer.

Attention

Damage to the bladder and tubing caused by disinfectants

- Before washing, carefully remove the bladder and tubing from the cuff sleeve.
- Avoid disinfecting the cuff bladder and connected rubber tubing.
- The bladder and tubing can be damaged by disinfectants.
- Wipe down the bladder with lukewarm water and add a mild detergent, if necessary.
- Ensure that no water enters the tube opening.

1. When cleaning the cuff sleeve, bladder and tubing, use only mild detergents in lukewarm water without fabric softener.
2. Read the safety instructions carefully and observe them closely during cleaning.

Maintenance plan

Attention

Damage to device

- Do not open the casing. Once the device is opened, all warranties will lapse.

Weekly Maintenance

Analysis review:

1. Review the print-out of your measurement analysis for:
 - Correctly entered times and intervals in accordance with the log.
 - Times of day/night transitions.
 - Correct standard values (nocturnal decrease).
2. Check the device, cuff and the cuff tubing for superficial soiling and clean it as specified in the **Cleaning** section.
3. Check the cuff and the cuff tubing for superficial damage. In the event of damages return it to your Welch Allyn specialist for inspection.

Checking battery voltage:

Always use fully charged or new batteries.

The battery voltage appears on the display of the ABPM 7100 for approximately 3 seconds after the device is switched on. The battery voltage must be at least 2.6 volt to ensure a 24-hour measurement.

Maintenance every 2 years

As proof of continuous compliance to “Basic Requirements” pursuant to Directive 93/42/EEC, the ABPM 7100 must be subjected to metrological checks every two years. In certain countries, this requirement may be regulated by national laws or regulations.

Welch Allyn offers to provide metrological checks and the servicing comprising of the following:

- Metrological monitoring.
- Software updates (if required)
- Functional check: Electronics, pump and pneumatic circuit.

Troubleshooting

Attention

Damage to device

- Do not open the casing. Once the device is opened, all warranties will lapse.

Basic error sources

The following may cause error measurements or unintended events:

- The patient’s arm movement during measurement
- Device deactivation (e.g. at night)
- Incorrect cuff size
- Cuff displacement while wearing it
- Omitted successful initial measurement by the doctor
- Medication not taken
- Wrong log set by the user
- Empty, incorrectly charged or outdated batteries
- Kinked or knotted cuff tubing

Transmission error

The ABPM 7100 reviews the transmitted data to prevent errors. If an error occurs, “**E004**” will be shown on the display.

Checklist

Please review the following checklist for any errors occurring during the operation of the ABPM 7100. Many errors have simple causes:

- Check to see that all cables are connected correctly.
- Check to see whether the ABPM 7100 and the computer are switched on.
- Check to see whether the batteries have sufficiently voltage.

Note Some errors are combined with a continuous alarm for safety reasons. The continuous alarm can be cancelled by pressing any button. If there is residual pressure inside the cuff, open the cuff immediately.

24 - Troubleshooting

Error codes

Error description of the ABPM 7100

Error symptom	Possible cause	Remedy
Time and date are not updated following a longer period without power supply from power packs or batteries.	The internal buffer battery is depleted.	Date and time can be reset after every power pack or battery replacement. Send the device to your Welch Allyn specialist.
Measurement data can no longer be called up/displayed.	An error occurred during patient data storage.	Delete the respective patient (menu bar) and recreate it.
The connection between the ABPM 7100 and the PC is faulty.	The incorrect COM interface is set.	Set the correct interface in the service programs.
	Cable plug or socket is defective.	Inspect the plug and the socket on the ABPM 7100. Ensure that the pins are straight to guarantee contact.
	The ABPM 7100 is not in transmitting mode (the displays shows the time).	Switch the ABPM 7100 off and then on again without removing the connecting cable.
No patient number.	The ABP Monitor is not initialized, i.e. the patient number was not transferred during the preparation for a 24-hour measurement.	The patient number can also be transmitted after the measurement. This does not influence the measurement data.
No measurements were conducted during the nocturnal phase.	The battery packs or batteries were prematurely depleted.	The power packs or batteries may be defective (please contact your Welch Allyn specialist).
	The patient has switched off the ABPM 7100.	Draw the patient's attention to the urgency of a complete 24-hour measurement.
The display does not show "co" or "bt".	You are not in transmitting mode.	Communication via cable: Switch the ABPM 7100 off and then on again without pulling the plug.
		Communication via BT: Press and hold the START button and then press the DAY/NIGHT button. Select "bt" using the START button.
No automatic measurements will be performed.	No manual measurements performed after application.	Valid manual measurement must always be performed after the device has been positioned.
	Incorrect log set.	Set log 1 or 2.
The measurement interval does not meet your expectations.	Incorrect log set.	The programmed log does not correspond with the set log in the ABPM 7100. Check the log manually on the device.
	No manual measurements performed after application.	Conduct manual measurement to activate the set log

Error symptom	Possible cause	Remedy
Err 1	The patient displays severe arrhythmia.	ABP Monitor not applicable.
	The arm was moved during measurement.	Keep the arm still during measurement.
	Insufficient valid pulse rate detected.	Place the cuff on your arm again.
Err 2	The arm was moved during measurement.	Keep the arm still during measurement.
	Cuff does not fit the arm snugly.	Check the seating of the cuff and that of the device.
Err 3	Blood pressure beyond the measurement range.	Permanent notification renders the ABP Monitor unsuitable for the patient.
	Strong arm movements.	Keep the arm still during measurement.
	Problems with the pneumatics.	If the error persists permanently, send the device to your Welch Allyn specialist.
Err 4	Data transmission cable incorrectly inserted in ABP Monitor.	Insert the cable into the ABP Monitor correctly.
	Pins in the plug of the data transmission cable are mechanically damaged.	Check the plug to see whether the pins on the inside are damaged. If they are, contact your Welch Allyn specialist.
	Measurement value was not correctly transmitted.	Restart the transmission.
Err 5 bAtt	Power pack or battery voltage too low.	Replace the power packs or batteries.
	Power packs or batteries are defective.	The power pack or battery voltage is correct but “ bAtt ” is displayed during cuff inflation. Replace the power packs.
	Battery contacts are corroded.	Clean the battery contacts with a cotton cloth and a little alcohol.
Err 6 + Possible continuous alarm until a button is pressed.	Build-up 25 fair.	Check the cuff for a build-up of air or a kink in the tubing. If the cuff tubing is kinked, straighten the tubing. Otherwise send the device in for inspection immediately.
	Blood pressure cuff incorrectly connected.	Connect the cuff to the device.
	Leaky points in the cuff or the cuff tubing.	If necessary, replace the cuff.
Err 7	The memory of the blood pressure measurement device is full. (a maximum of 300 measurements and events can be stored)	Delete the data in the ABP Monitor but ensure that the data was stored on your PC first.
Err 8	Measurement cancelled with a pressed button.	
Err 9 + Possible continuous alarm until a button is pressed.	Residual pressure inside the cuff	Wait for the cuff to deflate completely.
	Zero point comparison was unsuccessful.	Send the device to your specialist for inspection immediately or directly to your Welch Allyn specialist.

26 - Troubleshooting

Error symptom	Possible cause	Remedy
Err 10 + Continuous alarm until a button is pressed.	Severe error caused by accumulated pressure outside the measurement process.	Send the device to your specialist for inspection and repair immediately or directly to your Welch Allyn specialist.
	These error messages all show a severe error in the program code.	
The analysis unit does not react to data transmission but the display shows "co".	Data transmission cable not correctly inserted in the PC. (also refer to Err 4)	Check whether the 9-pin plug of the data transmission cable is securely seated in the device's interface socket. (also refer to Err 4)
The ABPM 7100 measures every two minutes.	Log 9 is set in the ABPM 7100.	Set log 1 or 2.
The desired log cannot be set with the button combination.	The last patient's measurement values are still contained in the memory.	Delete the data in the ABP Monitor but ensure that the data was stored first.
The ABP Monitor cannot be switched on.	The battery packs or batteries were incorrectly inserted.	Reinsert either power packs or batteries and ensure correct polarity.
	Power pack or battery voltage too low.	Replace the power packs or batteries.
	Defective display.	Send the device to your specialist for repair or directly to your Welch Allyn specialist.
An error occurs during the first measurement.	The cuff size is not suitable for the patient's arm circumference.	Measure the patient's arm circumference and compare this with the imprint on the cuff. You may require a different cuff size.

Communication error ABPM 7100 Bluetooth Interface

Error symptom	Possible cause	Remedy
Code 1	Bluetooth® interface of the ABPM 7100 was not started correctly. Possible hardware fault.	Send the device to your Welch Allyn specialist for inspection.
Code 2	Bluetooth® interface of the ABPM 7100 could not be configured correctly. (Communication problem between ABPM 7100 and the Bluetooth® module.)	Try again. If the error persists, send the device to your Welch Allyn specialist for inspection.
Code 3	The status of the Bluetooth® interface of the ABPM 7100 could not be determined. (Communication problem between ABPM 7100 and the Bluetooth® module.)	Try again. If the error persists, send the device to your Welch Allyn specialist for inspection.
Code 4	The Bluetooth® interface of the ABPM 7100 has not yet been paired with the analysis software.	Reconnect the device via Bluetooth®.
Code 5	The Bluetooth® interface of the ABPM 7100 could not connect to the Bluetooth dongle on the computer.	Try again. If the error persists, send the device to your Welch Allyn specialist for inspection.
Code 6	The measurement value memory of the ABPM 7100 contains unsent blood pressure values.	These will be sent once further measurements have been performed.
Code 7	The ABPM 7100 is paired with a cell phone or GSM modem, which is technically incapable of transmitting measurement values, is outside the network range or is incorrectly configured.	Try again. If the error persists, contact your Welch Allyn specialist.

Limited Warranty

Welch Allyn warrants the product to be free of defects in material and workmanship and to perform in accordance with manufacturer's specifications for the period of one year from the date of purchase from Welch Allyn or its authorized distributors or agents.

The warranty period shall start on the date of purchase. The date of purchase is: 1) the invoiced ship date if the device was purchased directly from Welch Allyn, 2) the date specified during product registration, 3) the date of purchase of the product from a Welch Allyn authorized distributor as documented from a receipt from said distributor.

This warranty does not cover damage caused by: 1) handling during shipping, 2) use or maintenance contrary to labeled instructions, 3) alteration or repair by anyone not authorized by Welch Allyn, and 4) accidents.

The product warranty is also subject to the following terms and limitations: Accessories are not covered by the warranty. Refer to the directions for use provided with individual accessories for warranty information.

Shipping cost to return a device to a Welch Allyn Service center is not included.

A service notification number must be obtained from Welch Allyn prior to returning any products or accessories to Welch Allyn's designated service centers for repair. To obtain a service notification number, contact Welch Allyn Technical Support.

THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. WELCH ALLYN'S OBLIGATION UNDER THIS WARRANTY IS LIMITED TO REPAIR OR REPLACEMENT OF PRODUCTS CONTAINING A DEFECT. WELCH ALLYN IS NOT RESPONSIBLE FOR ANY INDIRECT OR CONSEQUENTIAL DAMAGES RESULTING FROM A PRODUCT DEFECT COVERED BY THE WARRANTY.

Service Policy

All repairs on products under warranty must be performed by Welch Allyn or by a service provider authorized by Welch Allyn. Unauthorized repairs will void the warranty. In addition, whether or not covered under warranty, any product repair should be performed exclusively by Welch Allyn or by a service provider that has been authorized by Welch Allyn.

If the product fails to function properly - or if you need assistance, service, or spare parts - contact the nearest Welch Allyn Technical Support Center.

Before contacting Welch Allyn, try to duplicate the problem, and check all accessories to ensure that they are not causing the problem. When calling, please be prepared to provide:

- Product name, model number, and serial number of your product.
- Complete description of the problem.
- Complete name, address and phone number of your facility.
- For out-of-warranty repairs or spare parts orders, a purchase order (or credit card) number.
- For parts orders, the required spare or replacement part numbers.

If your product requires warranty, extended warranty, or non-warranty repair service, please call first the nearest Welch Allyn Technical Support Center. A representative will assist you troubleshooting the problem and will make every effort to solve it over the phone, avoiding potential unnecessary return of your product.

In case a return cannot be avoided, the representative will record all necessary information and will provide a Return Material Authorization (RMA) number, as well as the appropriate return address. An RMA number must be obtained prior to any return.

If you have to return your product for service, follow these recommended packing instructions:

- Remove all hoses, cables, sensors, power cords, and other accessories (as appropriate) before packing, unless you suspect they are associated with the problem.
- Wherever possible use the original shipping carton and packing materials.
- Include a packing list and the Welch Allyn Return Material Authorization (RMA) number.

It is recommended that all returned goods be insured. Claims for loss or damage to the product must be initiated by the sender.

EMC Guidelines and Manufacturer Declaration

Table 1 – Guidelines and Manufacturer Declaration
Electromagnetic emission for all ME devices and ME systems

Guidelines and Manufacturer Declaration - Electromagnetic emissions		
The ABPM 7100 is intended for operation in an electromagnetic environment as specified below. The customer or user of the ABPM 7100 should ensure its operation in such an environment.		
Emission measurement	Compliance	Electromagnetic Environment Guideline
HF Emissions according to CISPR 11	Group 1	The ABPM 7100 utilizes HF power for its internal function only. Its HF emission is therefore very low and it is improbable that neighbouring electronic device experience any interference.
HF Emissions according to CISPR 11	Class B	The ABPM 7100 is suitable for use in other facilities than the living area and those immediately connected to the public supply network, which also supplies buildings used for residential purposes.
Emission of harmonics according to IEC 61000-3-2	Not applicable	
Emission of voltage fluctuations/ flickers according to IEC 61000-3-3	Not applicable	

Table 2 – Recommended safety distances between portable and mobile HF communication devices and the ME device or ME system for ME devices or ME systems, which are not life-supporting


Recommended safety distances between portable and mobile HF telecommunication devices and the ABPM 7100			
The ABPM 7100 is intended for operation in an electromagnetic environment, in which disturbance variables are controlled. The customer or user of the ABPM 7100 can assist in preventing electromagnetic interferences by adhering to the minimum distance between portable and mobile HF telecommunication devices (transmitters) and the ABPM 7100 – depending on the output power of the device, as specified below.			
The transmitter's nominal capacity W	Safety distance, depending on transmission frequency m		
	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz
	$d = 1.2 \sqrt{P}$	$d = 1.2 \sqrt{P}$	$d = 2.3 \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23
For transmitters whose nominal capacity is not covered by the table above, the recommended safety distance d can be calculated in meters (m) using the formula given in the respective column, whereby P is the maximum nominal capacity of the transmitters in Watts (W) as per the information given by the transmitter's manufacturer.			
NOTE 1	At 80 MHz and 800 MHz, the higher frequency range applies.		
NOTE 2	These guidelines may not be applicable in all cases. The distribution of electromagnetic variables is influenced by the absorption and reflection of buildings, objects and people.		

**Table 3 – Guidelines and Manufacturer Declaration
Electromagnetic immunity – for all ME devices and ME systems**

Guidelines and Manufacturer Declaration - Electromagnetic immunity			
The ABPM 7100 is intended for operation in an electromagnetic environment as specified below. The customer or user of the ABPM 7100 should ensure its operation in such an environment.			
Immunity tests	IEC 60601-test levels	Compliance levels	Electromagnetic Environment - Guidelines
Electrostatic discharge (SD) according to IEC 61000-4-2	± 6 kV Contact discharge ± 8 kV Air discharge	± 6 kV Contact discharge ± 8 kV Air discharge	Floors should consist of wood or cement or ceramic tiles. If the floor consists of synthetic materials, relative humidity must be at least 30%.
Rapid transient electrical disturbance/bursts according to IEC 61000-4-4	± 2 kV for mains cables ± 1 kV for input and output lines	Not applicable ± 1 kV for input and output lines	The ABPM 7100 does not have an AC power supply.
Surges according to IEC 61000-4-5	± 1 kV Line-to-line voltage ± 2 kV Line-to-earth voltage	Not applicable Not applicable	The ABPM 7100 does not have an AC power supply.
Voltage drops, short interruptions and fluctuations in supply voltage according to IEC 61000-4-11	< 5% UT (> 95% drop in UT) for a 1/2 period 40% UT (60% drop in UT) for 5 periods 70% UT (30 % drop in UT) for 25 periods < 5% UT (> 95% drop in UT) for 5 seconds	Not applicable Not applicable Not applicable Not applicable	The ABPM 7100 does not have an AC power supply.
Magnetic field in supply frequency (50/60 Hz) according to IEC 61000-4-8	3 A/m	3 A/m	Magnet fields in network frequency should correspond with the typical values found in business and hospital environments.
NOTE UT is the AC voltage before the application of the test levels.			

30 - EMC Guidelines and Manufacturer Declaration

Table 4 – Guidelines and Manufacturer Declaration
Electromagnetic immunity for ME devices or ME systems that are not life-supporting

Guidelines and Manufacturer Declaration - Electromagnetic immunity			
The ABPM 7100 is intended for operation in an electromagnetic environment as specified below. The customer or user of the ABPM 7100 should ensure its operation in such an environment.			
Immunity tests	IEC 60601-test levels	Compliance level	Electromagnetic Environment - Guidelines
<p>Conducted disturbance variables according to IEC 61000-4-6</p> <p>Radiated disturbance variables according to IEC 61000-4-3</p>	<p>Effective value 3 V 150 kHz to 80 MHz</p> <p>3 V/m 80 MHz to 2.5 GHz</p>	<p>3 V</p> <p>3 V/m</p>	<p>Portable and mobile RF communications equipment should be used no closer to any part of the ABPM 7100, including cables, than the recommended safety distance calculated from the equation applicable to the frequency.</p> <p>Recommended safety distance:</p> $d = 1.2 \sqrt{P}$ <p>$d = 1.2 \sqrt{P}$ for 80 MHz to 800 MHz</p> <p>$d = 2.3 \sqrt{P}$ for 800 MHz to 2.5 GHz</p> <p>with P as the transmitter's nominal value in watt (W) according to the transmitter manufacturer's information and d as the recommended safety distance in meters (m).</p> <p>The field intensity of stationary radio transmitters, as determined by an electromagnetic on-site^a survey, should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol.</p> 
NOTE 1	At 80 MHz and 800 MHz, the higher frequency range applies.		
NOTE 2	These guidelines may not be applicable in all cases. The distribution of electromagnetic variables is influenced by the absorption and reflection of buildings, objects and people.		
<p>^a Theoretically, the field intensity of stationary transmitters such as base stations of cellular telephones and mobile landline devices, amateur radio stations, AM and FM radio and TV stations, cannot be accurately determined in advance. To assess the electromagnetic environment regarding the stationary transmitter, a study on the location's electromagnetic phenomena should be considered. If the measured field strength in the location in which the ABPM 7100 is used exceeds the applicable RF compliance level above, the ABPM 7100 should be monitored to verify normal function. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the ABPM 7100.</p> <p>^b Over the frequency range of 150 kHz to 80 MHz, field intensity should be less than 3 V/m.</p>			

Patient Information - operation of the ABPM 7100

Safety instructions

Warning

Risk of strangulation posed by the shoulder strap and cuff tubing.

- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Do not place the shoulder strap and cuff tubing around the patient's neck.
- Always place the cuff tubing under the outer clothing (even at night).
- Switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.
- Measurement can be interrupted at any stage by pushing a random button. This deflates the cuff and the device can be removed.

Warning

Poor circulation caused by continuous cuff pressure.

- Do not kink the connecting tubing.
- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Ensure the correct placement of the shoulder strap and cuff tubing.
- Always place the cuff tubing under the outer clothing (even at night).
- Switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Warning

Placement and inflation of the cuff over a wound may lead to further injuries.

Placement and inflation of the cuff on any limb with an intravascular access or under intravascular treatment or an arteriovenous (A-V) shunt may lead to temporary interruption of circulation and therefore to further patient injury.

Placement and inflation of the cuff on the arm at the side of a breast amputation may lead to further injury.

- Switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Warning

If the patient is wearing an additional ME device on the same limb for monitoring purposes, the placement and inflation of the cuff may trigger the temporary loss of the existing ME device's function.

The operation and use of the automated non-invasive blood pressure monitoring device may lead to longer impaired blood circulation in the patient or respective limb.

- Switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

Warning

Poor circulation due to overly frequent measurements.

- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Switch off and remove the device immediately and to inform the doctor in the event of experiencing any pain.

32 - Patient Information - operation of the ABPM 7100

Caution

Risk of injury caused by incorrect application of the cuff.

- If the patient has limited cognitive abilities, the device may only be used under supervision.
- Ensure that neither the shoulder strap nor the cuff tubing can ever wrap around the patient's neck. Always place the cuff tubing under the outer clothing (even at night).
- Place the device in such a way that, while the cuff is inflated, the tubing is not compressed or kinked, especially during sleep.
- Petechiae, haemorrhages or subcutaneous haematoma may occur in some patients.
- Switch off and remove the device immediately and inform the doctor in the event of experiencing any pain.

Attention

Damage to device

- Do not open the casing. Once the device is opened, all warranties will lapse.

Attention

Damage to device

- Do not wear the ABPM 7100 while showering. If you suspect that liquid has entered the device while cleaning or using it, the device shall no longer be used on the patient.
- If the device was exposed to moisture, switch off the device and remove the batteries.
- The device may not be operated around MRI scanners or in the immediate vicinity of other medical electrical equipment.
- During a defibrillator discharge, the device shall not be in contact with the patient. Such a discharge may damage the ABPM 7100 and cause it to display incorrect values.
- Measurement can be interrupted at any stage by pushing a random button. This deflates the cuff and the device can be removed.

Attention

Measurement errors

- Although the ABPM 7100 fulfills all EMC standard requirements, it should not be exposed to any strong electromagnetic fields, as this may lead to malfunctions outside the limit values.
- The cuff tubing between the ABPM 7100 and the cuff may not be knotted, compressed or pulled apart.
- The cuff connection must always engage with an audible "*CLICK*". A loose connection between the tubing and the device leads to measurement errors.

Note

- Severe malfunctions are indicated by a continuous beep.
- In the event of a continuous beep, switch off the device, remove the cuff and inform your doctor.

24-hour measurement

1. Before a 24-hour measurement, go through these instructions together with your doctor.
2. Let your doctor explain possible hazards in detail on the basis of the warnings above.
3. Ensure that you have understood all functions and observable points.



Safety:

For your own safety during the following steps, please observe the safety instructions at the start of this chapter.

Position of the cuff

It is imperative that the artery symbol is positioned on the brachial artery. If the cuff is aligned correctly, the metal bar will lie on the outside of the upper arm (on the elbow side), whereby the fabric bag must cover the skin under the metal bar.

The Buttons



ON/OFF

The **ON/OFF** button turns on and off the ABPM 7100 when the button is pressed for more than 2 seconds.



START

The **START** button serves to

- initiate the automatic protocol.
- trigger a measurement in addition to the automatic protocol.



DAY/NIGHT

The **DAY/NIGHT** button is used to differentiate between waking and sleeping phases during the measurement. Press the **DAY/NIGHT** button immediately before going to bed and upon waking.



EVENT

Press the **EVENT** button to record an event which may affect the blood pressure and to trigger an additional measurement. Note the reason for pressing the **EVENT** button in the event log.

Measurement process

During the first measurement, the cuff is inflated in increments, to determine the cuff pressure required to measure the systolic blood pressure value. This maximum required inflation pressure is stored and applied by direct inflation during the subsequent automatic measurements. The patient should stay calm during the measurement process, until the measurement is completed. Allow your arm to hang loose, or place your lower arm loosely on the table or on a support whilst sitting. Avoid any movement! In the event of a failed measurement a new measurement is performed automatically according to the measurement process described above.

Cancelling a measurement

A measurement will be cancelled by pressing any buttons during the measurement process causing the cuff to be quickly deflated automatically. The LCD display will then show “-**StoP**-” and the ABPM 7100 will beep 5 times. This cancellation will be stored in the measurement value table under **Cancel**.

Reorder No. 106169
Material No. 722592
DIR 80019691 Ver. A

