

**deSoutter**  
MEDICAL

Total Power, Versatility and Reliability



**orthodrive®**

[www.de-soutter.com](http://www.de-soutter.com)



*High performance motor offers unrivalled large bone reaming and sawing capability*

The **orthodrive®** system is the latest innovation in high performance large bone power tools from De Soutter Medical. The MBQ high torque rotary handpiece provides the surgeon with ample power even for the most demanding acetabular reaming procedures. The new drill/ream handpiece produces 20% more reaming torque and the powerful dedicated sagittal saw cuts through the hardest of bone with ease and precision.

## POWER

- Through continual engineering refinements, including our fourth generation hybrid motor controller and intelligent battery charging technology, we have increased power by a further 20%
- Li-Ion battery technology assures class leading power to weight ratio.

## VERSATILITY

- Exceptional system versatility. 15 different rotary attachments, plus 8 different function attachments, including: reciprocating and mini sagittal saws, wire/pin driver and a unique pulse lavage system.
- In-built, quickchange drill/ream gearbox eliminates dedicated reaming attachments, reducing cost and weight.
- Optional twin trigger rotary handpiece with instant forward, reverse and unique oscillating drill/tap modes.

## RELIABILITY

- Proven system reliability enhanced through continuous engineering developments, including a hermetically sealed hybrid motor controller and unique intelligent battery charging technology.
- New high performance surface finish protects against the most aggressive detergents whilst offering a highly scratch resistant barrier.

Whether your preference is for battery - either aseptic or sterile - or indeed corded power, the **orthodrive®** system provides total OR flexibility.









**MBQ-701 Double trigger handpiece fitted with reamer attachment**

**MBQ-700 Single trigger handpiece fitted with A.O. drill chuck**



*See pages 6 and 7 for drilling, reaming, sawing, wire driving and pulse lavage attachment options.*

	Model	Order code	Drilling Speed	Reaming Speed	Maximum reaming torque	Cannulation
	MBQ-700	1289774	0 - 1350 rpm	0 - 300 rpm	20 Nm	Ø 4.4 mm
	MBQ-701	1289884	0 - 1350 rpm	0 - 300 rpm	20 Nm	Ø 4.4 mm



**DBK-701**  
Sagittal saw  
shown with aseptic  
Li-Ion housing and  
battery

**DBC-700**  
Reciprocating saw  
shown with sterile  
Li-Ion battery

*Ask for Stericut®  
catalogue  
for range of sterile  
packed saw blades.*

	Saw Type	Model	Order code	High speed
	Reciprocating	DBC-700	1290444	0 - 13500 cpm
	Sagittal (S89)	DBK-701	1290004	0 - 13500 cpm
	Sagittal (S87)	DBK-700	1289994	0 - 13500 cpm



**P30-700**

*Pump kit attached to MBQ-700 handpiece with a simple "click in" design*

The ECO Pulse doesn't require batteries, motors or electronics and therefore typically saves up to **55%** in clinical waste, per procedure.

A smaller footprint allows less packaging to be used, this results in a **reduction in valuable storage space by 2.5 times.**

**PUMP KITS**

Order code	Type	Supplied nozzle
P30-700	General purpose with suction	P30-121
P30-701	Hip arthroplasty with suction	P30-121/P30-123
P30-702	Hip arthroplasty brush with suction	P30-124/P30-121
P30-704	Hip arthroplasty without suction	P30-120/P30-122
P30-705	Bone cleaning (high pressure) with suction	P30-125
P30-707	Hip arthroplasty (high pressure) with suction	P30-125/P30-123
P30-708	Hip arthroplasty (high pressure) without suction	P30-125/P30-122
P30-709	Bone cleaning (fan spray) with suction	P30-127
P30-710	Hip arthroplasty (fan spray) with suction	P30-123/P30-127

**NOZZLES**

	Order code	Type	Useable length
	P30-125	Bone cleaning - high pressure with suction	6cm
	P30-127	Bone cleaning - fan spray with suction	6cm
	P30-121	Shower spray with suction	6cm
	P30-120	Soft tissue without suction	10cm
	P30-123	Femoral canal with suction	25cm
	P30-124	Femoral brush with suction	25cm
	P30-122	Femoral canal without suction	25cm

Flow rate Variable 0 - 1.0 litres/min


Pulse rate Variable 0 - 1350 pulses/min


\* Based on AB-703 Performance figures dependent on battery type.







## MBQ DRIVES AND ATTACHMENTS

	Order code	Model	Chuck description	Final drive speed	Cannulation Ø mm
DRILLING AND REAMING ADAPTORS					


	15360	-	0.5 - 4.0 mm (5/32") keyed	As handpiece	3.1
	15370	-	0.5 - 7.4 mm (1/4") keyed	As handpiece	4.2


	15380	-	0.3 - 3.0 mm (1/8") keyless	As handpiece	3.0
	15390	-	0.3 - 7.4 mm (1/4") keyless	As handpiece	4.2

	15480	-	Small A.O. Synthes/ASIF	As handpiece	2.5
---	-------	---	-------------------------	--------------	-----

	15490	-	Large A.O. Synthes/Protec	As handpiece	4.2
	15530	-	Hudson	As handpiece	4.2
	15540	-	Trinkle	As handpiece	4.2
	15460	-	Zimmer	As handpiece	4.2
	15470	-	Hudson/Zimmer combi	As handpiece	4.2
	17130	-	Harris	As handpiece	4.2
	19330	-	DHS/DCS	As handpiece	4.2
	18190	-	Stryker Triathlon®	As handpiece	2.3
	17930	-	Aesculap Taragon TX	As handpiece	4.2

### SAWING ATTACHMENTS


	15560	NQ-700	Mini sagittal saw Uses S88 type blades	0 - 24500 cpm	-
--	-------	--------	---	---------------	---

	18280	CQ-710	Reciprocating saw. Uses S84 type blades	0 - 12300 cpm	-
	18390	CQ-711	Reciprocating saw. Uses S22 type blades	0 - 12300 cpm	-


### RADIOLUCENT DRILL ADAPTOR

	15510	TQ-700	Accepts radiolucent drill attachment	As handpiece	-
---	-------	--------	--------------------------------------	--------------	---

### 90° RADIOLUCENT DRILL ATTACHMENT

	17710		Uses D56 type drill bits	As handpiece	-
---	-------	--	--------------------------	--------------	---

### WIRE/PIN DRIVERS

	17240	WQ-710	Wire/Pin Driver 0.7 - 4.0 mm cap.	As handpiece	4.0
	17940	WQ-701	Variable Wire Driver 0.6 - 2.0 mm cap.	As handpiece	2.0
	17860	PQ-701	Variable Pin Driver 1.8 - 3.2 mm cap.	As handpiece	3.2
	17950	PQ-702	Variable Pin Driver 2.6 - 4.0 mm cap.	As handpiece	4.0

Ask for Stericut® catalogue for range of sterile packed saw blades, drill bits, rotary burs and K wires.



**BC-700 Battery charger**  
shown with aseptic batteries.

The BC-700 Battery Charger requires minimal user input, simply insert battery and the BC-700 will charge, recover and recondition it to its optimum state.



**The large touch-screen display** with intuitive menu navigation, displays instant battery charge status and condition. Includes a multi-lingual display and other user settings.

The BC-700 also offers a **unique healthcheck facility**, profiling the condition of an individual battery, automatically reconditioning it, to provide optimal capacity. Alternatively, a manual healthcheck can be easily selected.

The BC-700 Battery Charger accepts all De Soutter Medical batteries, ensuring total backwards compatibility.



BC- 700 MULTI BAY CHARGER AND MODULE OPTIONS

Charger Order code	Charger module description	Charger module Order code	For battery type
15700	Without charger module fitted	-	-
16230	With 4-off CA-703 charger modules	18120	 AB-703/AB-704
16220	With 4-off CA-700 charger modules	15750	 SB-703/SB-704

Mains power cordsets - 3 metres long. Please specify	-	Order code
UK		8290
Euro		8300
US/Canada		8310
Australia/New Zealand		8320
Swiss		14700

BC-700 charger modules can be quickly and easily changed to suit the batteries being charged.



SINGLE BAY CHARGERS

Model	Battery type Order code	For battery type
BC-706	17990	SB-703/SB-704
BC-707	18000	AB-703
BC-708	18010	AB-704



# LITHIUM ION ASEPTIC BATTERIES

## Why choose Lithium Ion battery technology

De Soutter Medical Lithium Ion batteries offer the surgical team and sterile services a number of tangible benefits over conventional battery technology.

Key advantages are:

- Li-Ion technology improves battery 'capacity to weight' ratio by around 25% when compared to equivalent NiMH or NiCad chemistries.
- Significantly extends the useable run time of the power tool.
- Eliminates 'memory effect' normally associated with other battery types.
- De Soutter Medical sterile Li-Ion technology readily withstands repeated machine washing and autoclaving.
- Unique **fuel gauge** incorporated within De Soutter Medical Li-Ion batteries allows the surgical team to instantly assess the battery capacity - both before and during use.

## ASEPTIC LITHIUM ION BATTERIES

### SMALL BATTERY



Battery pack



Aseptic housing



Aseptic shield

Model	Description	Order code	Voltage	Capacity mAh
AB-703	Battery pack	17910	10.95V	2000
AH-703	Battery housing	18080	-	-
AS-703	Battery shield	18100	-	-

### LARGE BATTERY



Battery pack



Aseptic housing



Aseptic shield

Model	Description	Order code	Voltage	Capacity mAh
AB-704	Battery pack	17920	14,6V	2000
AH-704	Battery housing	18090	-	-
AS-704	Battery shield	18110	-	-

STERILE LITHIUM ION BATTERIES

SMALL BATTERY



Model	Order code	Voltage	Capacity mAh
SB-703	17210	13,2V	1100

LARGE BATTERY



Model	Order code	Voltage	Capacity mAh
SB-704	17220	13,2V	2500



CORDED DIRECT POWER



**PS-701 Power supply**  
Order code 17870

Connected to MBQ-701  
via corded module.

New slim line design.  
Requires no set-up. Simply connect  
to the CM-700 Corded Module for  
instant use.

Mains power cordsets  
3 metres long.  
Please specify from options below.

UK	Order code 8290
Euro	Order code 8300
US/Canada	Order code 8310
Australia	Order code 8320
New Zealand	Order code 8320
Swiss	Order code 14700

**CM-700  
Corded module**  
Order code  
15720




# STERILISATION CONTAINERS

*Aseptic tray*



*Sterile tray*

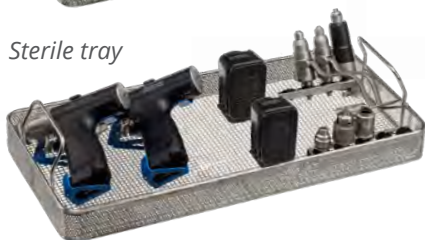

*Half size  
battery  
insert tray*

## ASEPTIC BATTERY STERILISATION SYSTEM

Set up - optional tray contents	Standard wire container and insert tray	Standard insert tray
3 off handpieces 9 off attachments 2 off aseptic battery housings and shields	Order code 18150	Order code 17960
2 off handpieces 9 off attachments 2 off aseptic battery housings and shields	Order code 18160	Order code 17970

## STERILE BATTERY STERILISATION SYSTEM

Set up - optional tray contents	Standard wire container and insert tray	Standard insert tray
3 off handpieces 9 off attachments 2 off sterile batteries	Order code 16250	Order code 17110
2 off handpieces 9 off attachments 2 off sterile batteries	Order code 19140	Order code 19040

## WIRE CONTAINERS

Standard wire container  
Order code 19010  
Accepts standard insert tray

Half size wire container  
Order code 19000  
Accepts half size insert tray

Set up - optional tray contents	Half size wire container and insert tray	Half size Insert tray
2 off SB-70x batteries	Order code 17390	Order code 17400

## STERILE FILTERED CONTAINERS

Accepts standard  
insert tray.



Standard sterile filtered container  
Order code 19250

PTFE sterile filtered container  
Order code 19810

Lid gasket - Order code 670023

Accepts half size  
insert tray.



Half size sterile  
filtered container  
Order code 19370  
Lid gasket - Order code 669903

Filter paper  
single use  
pack of 100  
Order code 657243



Security tag  
pack of 100  
Order code  
657253



Sterile indicator label  
pack of 100  
Order code  
657263



PTFE filter cassette  
Order code 671073



**UNITED KINGDOM (INTERNATIONAL HEAD OFFICE)**

De Soutter Medical Limited  
Halton Brook Business Park, Weston Road,  
Aston Clinton, Aylesbury, Bucks HP22 5WF  
Tel: +44 (0) 1296 634 000  
Fax: +44 (0) 1296 634 033  
Email: [info@de-soutter.com](mailto:info@de-soutter.com)

**AUSTRALIA**

De Soutter Medical Australia Pty Ltd  
2/12-14 Apollo Drive, Hallam  
Victoria 3803  
Tel: +61 (0) 3 9702 4441  
Fax: +61 (0) 3 9702 4484  
Email: [australia@de-soutter.com](mailto:australia@de-soutter.com)

**AUSTRIA**

De Soutter Medical Austria  
Zweigniederlassung Österreich  
Dietrichsteingasse 10, A-3400 Klosterneuburg  
Tel: +43 (0) 676 96 71 770  
Fax: +43 (0) 2243 21 656  
Email: [austria@de-soutter.com](mailto:austria@de-soutter.com)

**BELGIUM**

De Soutter Medical Belgium  
Bessemmerstraat 14  
3620 Lanaken  
Tel: +32 (0) 89/47 15 37  
Fax: +32 (0) 89/70 12 19  
Email: [belgium@de-soutter.com](mailto:belgium@de-soutter.com)

**FRANCE**

De Soutter Medical France  
1 252 Avenue du Parc des Expositions  
33260 La Teste de Buch  
Tel: +33 (0) 5 56 54 89 36  
Fax: +33 (0) 9 70 61 37 60  
Email: [france@de-soutter.com](mailto:france@de-soutter.com)

**GERMANY**

De Soutter Medical Deutschland  
Bahnhofstraße 4  
D-66625 Nohfelden  
Tel: +49 (0) 68 52-99 12 46  
Fax: +49 (0) 68 52-99 12 47  
Email: [deutschland@de-soutter.com](mailto:deutschland@de-soutter.com)

**ITALY**

De Soutter Medical Italia  
Località Fornace SNC  
27022 Casorate Primo - PV  
Tel: +39 (0) 2 9009 4098  
Fax: +39 (0) 2 9009 2673  
Email: [italy@de-soutter.com](mailto:italy@de-soutter.com)

**NETHERLANDS**

De Soutter Medical Nederland  
Gelderlandhaven 2X  
3433 PG Nieuwegein  
Tel: +31 (0) 85 0491480  
Fax: +31 (0) 85 0491489  
Email: [nederland@de-soutter.com](mailto:nederland@de-soutter.com)

**UNITED STATES OF AMERICA**

De Soutter Medical USA, Inc.  
224 Rolling Hill Road, Suite 12A  
Mooresville, NC 28117  
Tel: +1 (704) 655 9040  
Fax: +1 (704) 987 2035  
Email: [usa@de-soutter.com](mailto:usa@de-soutter.com)

**deSoutter**  
MEDICAL

For further information on de Soutter Medical powered instruments  
and cast saw systems please contact the address in your area